

# Signia Pocket Guide

May 2024 – Unleash the Power of Conversation  
with Signia Integrated Xperience

[veterans.signiausa.com](https://veterans.signiausa.com)



Pure Charge&Go T IX



Styletto IX and  
CROS Styletto IX



Inzio IX IIC/CIC



signia

# Table of Contents

## Technology

Signia Integrated Xperience .....	1
New Study .....	2
Signia Signature Features .....	3
Signia Signature Audiological Features .....	4
Be Active, Be Healthy, Be Brilliant with My WellBeing. ....	6
Signia App for Patients.....	7
Signia TeleCare – 3.0.....	8
Connex 9.11.10 UPDATE.....	9

## RIC

Styletto IX .....	11
Pure Charge&Go T IX .....	12
Pure 312 AX .....	13
Styletto AX .....	14
Signia RIC products and fitting ranges.....	15
miniReceiver 3.0.....	16
EarWear 3.0 .....	17
RIC Earmolds.....	18

## Custom

Insio 7IX CIC/IIC.....	19
Insio 7IX CIC/IIC Fitting Ranges .....	20
Insio Charge&Go AX Customs.....	21
Insio Charge&Go AX Customs Fitting Ranges .....	22
Insio Nx Custom .....	23
Insio Nx Custom Fitting Ranges .....	24

## BTE

Motion Charge&Go X .....	25
Motion Charge&Go X Fitting Ranges.....	26
Motion 13 and 13 P Nx.....	27
BTE Earmolds .....	28
Thin Tube 3.0 .....	29

## CROS & Accessories

CROS.....	31
CROS — RICs .....	32
StreamLine™ Mic .....	33
StreamLine™ TV.....	36
miniPocket™ .....	39

## Chargers

Charger Compatibility .....	41
Rechargeability .....	44
Portable IX/AX Charger .....	45
Insio Charge&Go Charger .....	46
Inductive Charger II .....	47
Styletto IX/AX Charger .....	48
Dry&Clean Charger.....	49
How to order a Charger in ROES.....	51

## Signia Quick Reference

Signia Nomenclature.....	53
Where to Find Serial Numbers .....	54
How to Measure for Signia 3.0 Receivers .....	55
Signia Earwear 3.0 Parts.....	56
Signia Wax Guard Guide .....	57
Earmold Ordering Tips .....	58
Accessories Tips .....	59
Signia Fitting Guide.....	60
Tinnitus Fitting Guide .....	63
CROS Quick Fit Guide.....	64
CROS/BICROS Listening Check .....	65
How to Pair .....	66

ROES Codes .....	67
------------------	----

Technology



Signia Integrated Xperience

## Finally, a hearing aid designed for real conversations

**Meet the world's first hearing aid with a multi-stream architecture capable of accurately pinpointing multiple moving speakers in real-time conversation.**

Improved hearing is about more than just sound amplification. It is about empowering individuals to immerse themselves in their surroundings, create meaningful connections, and cherish every moment of life's experiences.

Leveraging our profound understanding of dynamic conversations, including the head movements of the wearer and changing positions of speakers, we crafted Signia Integrated Xperience. This groundbreaking technology ensures uninterrupted conversation engagement in noise, addressing a critical challenge many hearing aid wearers face.

Our commitment is to support wearers in leading active lives, unshackled from the limitations of traditional hearing aids and free to contribute and foster deeper connections.

**Experience the transformative power of real-time, dynamic conversation enhancement with Signia Integrated Xperience.**



New Study



# Industry-Leading Group Conversation Enhancement

**Signia Integrated Xperience boasts an impressive 4.1 dB SNR improvement over the nearest competitor in a noisy group conversation scenario.**

Our latest study spotlights RealTime Conversation Enhancement's unique ability to target and enhance multiple speakers, ensuring unparalleled access to conversation in noisy group settings.

## An Innovative Methodology: Providing Evidence Relevant to Real Conversations

We simulated a real group conversation between two alternating speakers in noise. We then measured the output SNR in the ear of a KEMAR manikin, comparing Signia IX to premium devices for four other prominent manufacturers. This revealed a remarkable 4.1 dB SNR improvement<sup>1</sup> over the nearest competitor—a testament to RealTime Conversation Enhancement's industry-leading ability to dynamically enhance multiple speakers in noise.

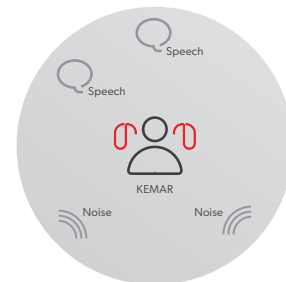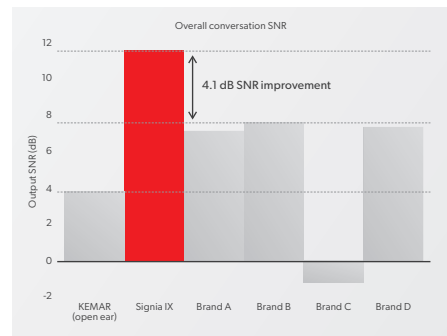

<sup>1</sup>Jensen NS, Wilson C, Kamkar Parsi H & Taylor B. 2023. Improving the signal-to-noise ratio in group conversations with Signia Integrated Xperience and RealTime Conversation Enhancement. Signia White Paper.

## Signia Signature Features



### **Auto EchoShield: Automatically better speech clarity and comfort even in reverberant situations**

Automatically analyzes the listening environment and adjusts to overcome problematic reverberation in every kind of room. Wearers will enjoy better speech clarity, less listening effort, better room perception, higher overall satisfaction, higher listening comfort, and better perceived speech intelligibility.

Visit [signia-pro.com/en-us/claims-data](https://signia-pro.com/en-us/claims-data) for more details



### **Own Voice Processing (OVP) 2.0: Unique solution to increase own voice comfort**

Can separate the wearer's own voice from other speech and noise for a natural sound, increasing first-fit acceptance to 88% and keeping patients engaged in conversations.

\*Not available with Inzio Charge&Go AX.



### **HandsFree with CallControl (iOS): Turn hearing aids into a headset**

The HandsFree with CallControl (iOS) features let the hearing aid wearer use their Signia IX and AX hearing aids for hands-free calls with their iPhone® in different environments — even in noisy ones. Wearers will enjoy answering the phone in any situation and never being limited in their active, busy lifestyle.

## Signia Signature Audiological Features

### Enhanced Audiological Features



#### **e2e Wireless 4.0: A great leap for true binaural hearing**

e2e 4.0 significantly boosts the speed of data transmission and a 60% reduction in latency provides inaudible shifts between environments. The synchronization of more binaural audio data than ever before ensures clearer speech processing with high-precision directional awareness.



#### **TruEar 360: Mimics the natural pinna effect**

The powerful processing engine allows for an exceptional TruEar 360 mic mode that reproduces and emphasizes the capabilities of the natural pinna-effect better than ever. TruEar 360 highlights the difference between front and back sources of sound to make localization of sound sources easier.



#### **Automatic Situation Detection : automatically responds to where and what you are doing**

Provides detailed information about the wearer's movement and environment. Analyzing and adapting to the acoustic and motion changes, these sensors deliver a clear understanding of speech from any direction and maximizing audibility as you move throughout the day.



## Signia Signature Audiological Features

### Enhanced Audiological Features



#### Enhanced feedback cancellation: No distractions for a smooth sound impression

e2e 4.0 also synchronizes feedback cancellation between the ears in the same phase, removing more artifacts for a smoother overall sound experience.



#### Direct Streaming: hassle-free listening to audio input

Signia Integrated Xperience and Augmented Xperience offer direct streaming for both iOS and ASHA-compatible Android™ smartphones. They handle the signal of a streaming input independently from other signals, ensuring that the sound is even clearer and easier to optimize.



#### Increased input dynamic range input: appreciated by music lovers

The microphones and the state-of-the-art analog input stage can handle signals up to 117dB with great sound quality. That's very close to the full capacity of human ears and ensures a clear, crisp signal, for scenarios like concerts.

## Be Active, Be Healthy, Be Brilliant with My WellBeing

With the new My WellBeing features in the Signia app, you can now see important health indicators empowering you to better maintain your physical activity.

Research shows that health trackers like My WellBeing help you stay active,\* which boosts your energy levels and confidence.\*\*

Get started today with Signia!



### My Steps

Staying physically active is shown to improve energy levels and confidence – and even slow down age-related hearing loss.\*\* The motion sensor in Signia AX hearing aids not only helps you hear better on the move, but also enables you to see exactly how many steps you take each day.



### My Activity

The use of a fitness tracker can help increase physical activity and keep you active.\* Signia My WellBeing is smart enough to know how much time you spend exercising — and together with your steps, provides an overview of your physical activity level.



### My WearTime

You can only take advantage of the many benefits your hearing aids provide when you are wearing them. My WellBeing monitors how well you are using your hearing aids based on how often you use them – helping you better understand your usage habits.



### My Conversations

Your social interactions play an important role in maintaining good hearing health. Signia's unique Own Voice Processing technology measures how much you are socially interacting with others, so you can know the duration and dynamics of social interactions.

My Conversations only available with OVP 2.0

\* Effectiveness of physical activity monitors in adults: systematic review and meta-analysis. BMJ, 376, e068047.

\*\* Fitness and Better Hearing, Part 1. Hearing Review, 21(3): 18-24



## Signia App for Patients

Beyond streaming and remote control capabilities, the Signia app is customized to display specific features based on the patient's hearing aid model and gives them freedom and support with access to:




**TeleCare™** for remote care support with you! Patients fitted with Signia TeleCare have a 19% higher acceptance rate than patients fitted without it.<sup>1</sup>



**Signia Assistant** for on-demand AI support. 93% stated that using the Signia Assistant improved their satisfaction for listening in difficult situations.<sup>2</sup>



**Signia Face Mask Mode** for an immediate boost in speech with the simple tap of a button. Facemasks have been shown to reduce audibility by up to 30%, and make lip reading impossible.<sup>3</sup>

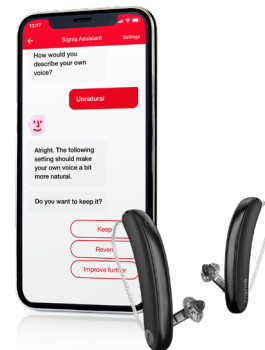

<sup>1</sup> Article, Matthias Froehlich, PhD, Eric Branda, AuD PhD, Dirk Apel, BS, 2020: User Engagement with Signia TeleCare: A Way to Facilitate Hearing Aid Acceptance

<sup>2</sup> Hoydal EH, Fisher R-L, Wolf V, Branda E, Aubreville M. Empowering the wearer: AI-based Signia Assistant allows individualized hearing care. Hearing Review. 2020;27(7):22-26.

<sup>3</sup> Goldin A, Weinstein BE, Shiman N. How do medical masks degrade speech perception? Hearing Review. 2020;27(5):8-9.

## Signia TeleCare – 3.0

Proven remote care support to remain virtually connected to your patients

Remain connected with patients — Communicate with patients through Signia’s safe and secure connection on any laptop, smartphone, or computer.


Higher patient satisfaction — More flexibility to support patient needs.

- In-person or remote activation of TeleCare for new and existing patients.
- Remote adjustments via Connexx.
- Compatible with multiple Signia platforms (Nx, X, AX, IX).
- Adjustments are sent to the patient in real-time (synchronously).



Patients fitted with Signia TeleCare have a 19% higher acceptance rate than patients fitted without it.\*

TeleCare Code  
no longer needed for  
easier Telecare Use

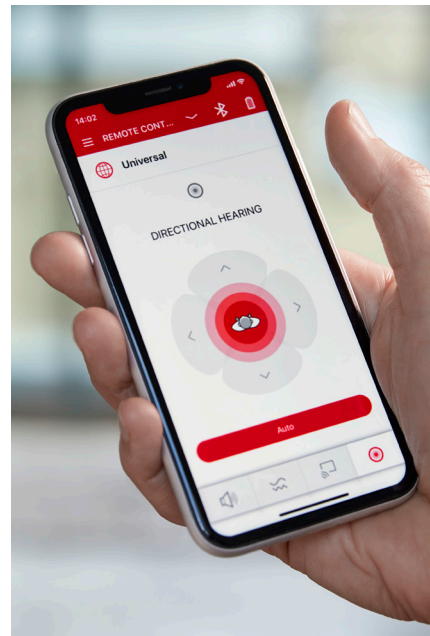

\* Article, Matthias Froehlich, PhD, Eric Branda, AuD, PhD, Dirk Apel, BS, 2020: User Engagement with Signia TeleCare: A Way to Facilitate Hearing Aid Acceptance

## Connexx 9.11.10 UPDATE

- Improved performance speed of software.
- Wireless hearing aid firmware updates are available via ConnexxAir wireless, only for non-Bluetooth® hearing aids. Firmware updates are only available when received from the notification center.



version **9.11.10** with SIFIT 9.11.10









# Styletto IX

The world's slimmest & most advanced SlimRIC



Styletto IX Charger



Styletto Dry&Clean Charger

20 hours on a single charge  
Portable Charge&Go Charger holds 5 charges

- Receivers: 3.0 S, M, P
- EarWear: EarWear 3.0

## Key features & benefits



RealTime Conversation Enhancement



Augmented Focus



Own Voice Processing 2.0



Signia Assistant



iPhone® & Android™ connectivity



Convenient rechargeability



Signia app



CROS compatibility



Sleek, stylish design



IP68



Black/Black Gloss



Black/Graphite



Black/Silver



Cosmic Blue/Rose Gold



Snow White/Rose Gold



Snow White/Silver



Snow White/Snow White Gloss



Beige/Beige Gloss



Deep Brown/Deep Brown Gloss



Forest Camo/Black



Desert Camo/Black

# Pure Charge&Go T IX

Unleash the Power of Conversation



IP68



Pure Charge&Go Charger



Portable Charge&Go Charger (Default)



Dry&Clean Charge&Go Charger

36 hours on a single charge  
Portable Charge&Go Charger holds 3 charges

- Receivers: 3.0 S, M, P, HP
- EarWear: EarWear 3.0



Black



Graphite



Dark Champagne



Silver



Pearl White



Fine Gold



Deep Brown



Sandy Brown



Rose Gold



Beige

## Key features & benefits



RealTime Conversation Enhancement



LE Audio Ready



Recharge on-the-go



CROS



Augmented Speech Understanding



Immersive Soundscape



Automatic Situation Detection



24/7 Signia Assistant



Auto EchoShield



Own Voice Processing 2.0



HandsFree with CallControl for iOS



My WellBeing



Integrated Telecoil



iPhone & Android Connectivity

## Pure 312 AX



IP68

- Receivers: miniReceivers 3.0 S, M, P, HP.
- Receivers are field replaceable
- EarWear: EarWear 3.0

### Key features & benefits



312 Battery



iPhone & Android Connectivity



24/7 Signia Assistant



CROS



Augmented Speech Understanding



Immersive Soundscape



Automatic Situation Detection



Auto EchoShield



Own Voice Processing 2.0



HandsFree with CallControl for iOS



My WellBeing



Black



Graphite



Dark Champagne



Silver Hero



Pearl White



Fine Gold



Deep Brown



Sandy Brown



Rose Gold



Beige

## Styletto AX

Exchangeable receivers, rechargeability, and *Bluetooth*® streaming in a one-of-a-kind design.

Available colors



IP68



Portable Charger (Default)  
Holds 4 days of charging.  
17 hours on a single charge.

Dry&Clean Charger

**Qi Wireless Technology<sup>1</sup>**  
Charger equipped with Qi Wireless Technology — enables wireless charging on any standard charging pad, including those built into cars and cafe counters.

### Key features & benefits



My WellBeing



iPhone & Android Connectivity



Recharge on-the-go



CROS



Augmented Speech Understanding



Immersive Soundscape



Automatic Situation Detection



24/7 Signia Assistant



Auto EchoShield



Own Voice Processing 2.0



HandsFree with CallControl for iOS

- Receivers: miniReceivers 3.0 S,M,P
- Earwear: Earwear 3.0



reddot design award  
winner 2019




<sup>1</sup> Wireless charging pad not included with purchase  
<sup>2</sup> Requires StreamLine Mic






# Signia RIC products and fitting ranges

## miniReceiver 3.0

### 6 Lengths

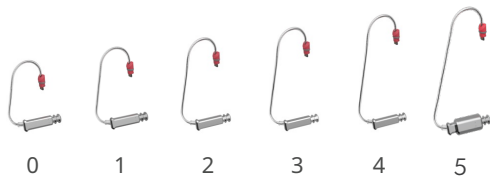

### 4 Power Levels

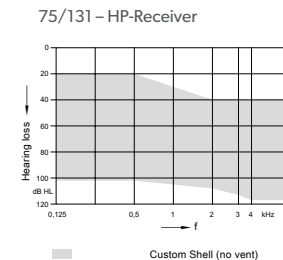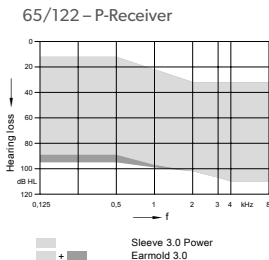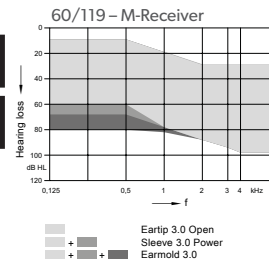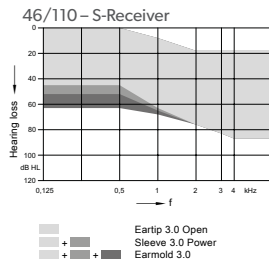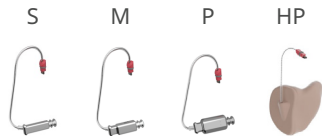

## miniReceiver 3.0

receiver types

Our new EarWear 3.0 offers a higher physical fit rate, improved cosmetics and better retention




miniReceiver 3.0 Cable Bending  
10 / 10 fit correctly – out of the box

Smaller plug for smaller RICs

Push-/pull system –  
no tool required

Smart Receiver Detection  
prevents misuse

6-Pin Connector for  
new features

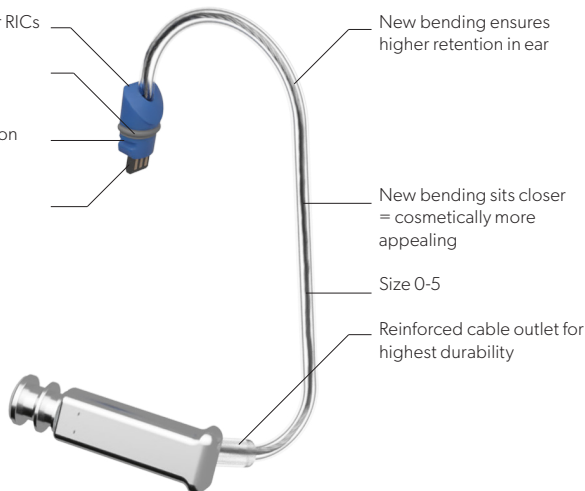

New bending ensures  
higher retention in ear

New bending sits closer  
= cosmetically more  
appealing

Size 0-5



Reinforced cable outlet for  
highest durability

## EarWear 3.0

This versatile new range is suitable for the vast majority of ear anatomies and hearing needs




Soft eartips allow for deep insertion with high comfort - also into small ear canals

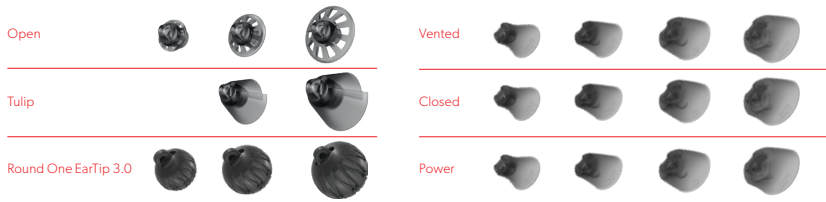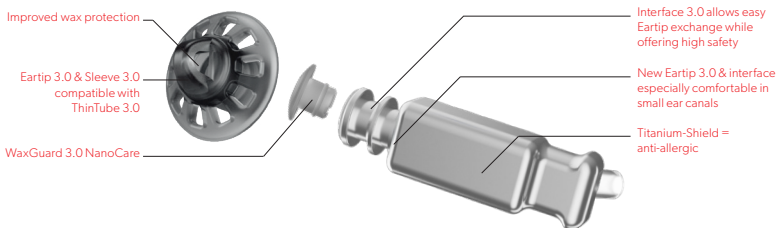- The new Power Sleeve 3.0 expands on the the acoustically reliable sleeve design, now offering a fully closed option.
- Secure, proven design makes Eartips 3.0 and Sleeves 3.0 much easier for wearers to exchange.
- All Earwear 3.0 items can be used together.



# RIC Earmolds

Bold text indicates default options

	RIC Earmold 3.0		RIC Embedded Earmold 3.0		RIC HP Earmold 3.0	
						
<b>Styles</b>	<b>RIC Hollow (Hard Material Only)</b> RIC Lock Hollow (Hard Material Only) Half Shell Hollow (Hard Material Only) Skeleton Hollow (Hard Material Only) RIC RIC Lock	Helix Lock Full Shell Half Shell Semi-Skeleton Skeleton	<b>RIC Solid</b> RIC Lock Solid Skeleton		<b>Canal</b> Canal Lock Helix Lock Skeleton Lock	
<b>Receivers</b>	Receiver 3.0 S (Size 0-5) Receiver 3.0 M (Size 0-5) Receiver 3.0 P (Size 0-5)		Receiver 3.0 S (Size 0-5) Receiver 3.0 M (Size 0-5) Receiver 3.0 P (Size 0-5)		Receiver 3.0 HP (Size 0-5)	
<b>Materials</b>	<b>Hard</b>	Soft	<b>Hard</b>		<b>Hard</b>	
<b>Colors</b>	Gloss and Matte Finishes Available <b>Clear</b> , Rose, Brown (Hard Material Only), Medium Brown, Dark Brown		Gloss and Matte Finishes Available <b>Clear</b> , Rose, Brown (Hard Material Only), Medium Brown, Dark Brown		<b>Beige</b> , Mocha, Brown, Red/Blue, Clear	
<b>Venting</b>	None Open (Hollow Only) Pressure (1.0MM) 1.2MM 1.6MM Large as Possible	SAV (2.0MM) SAV (2.5MM) SAV (3.0MM) Trench (hard material only) Semi-IROS	No Vent 1.0mm 1.2mm 1.6mm SAV (2.0MM)	SAV (2.5MM) SAV (3.0MM) As Large As Possible Trench Vent Semi-IROS	None Pressure (1.0 MM) 1.6 MM	
<b>Ear Canal Length</b>	Short, Medium, Long (Hard Material Only)		Short, Medium, <b>Long</b>		Short, Medium, Long, <b>Deep</b>	
<b>Indicators</b>	<b>Red/Blue Dots</b> Red/Blue Letters	Patient Name	<b>Red/Blue Dots</b> Red/Blue Letters	Patient Name	Large Dot Large R-L Letters	Patient Name
<b>Dexterity</b>	Removal Line		Removal Line		Removal Line	
<b>Wax Protection</b>	<b>NanoCare 3.0</b>		<b>NanoCare 3.0</b>		<b>Cerustop 3.0 XL</b>	



Custom



## Insio 7IX CIC/IIC

Engineered for conversations  
Filter out the noise  
Keep the connection  
Only with Signia Insio IX



### IIC:

- Faceplate:
- Black
  - Dark Brown
  - Mocha

- Shell:
- Blue
  - Red
  - Black
  - Dark Brown
  - Tan

### CIC:

- Faceplate:
- Black
  - Brown
  - Mocha

- Shell:
- Mocha
  - Brown
  - Black
  - Blue
  - Red

### Key features & benefits



Binaural OneMic directionality



e2ewireless



Dynamic Soundscape Processing 3.0



Convenient rechargeability






Signia app




IP68

# Insio 7IX CIC/IIC Fitting Ranges

## CIC

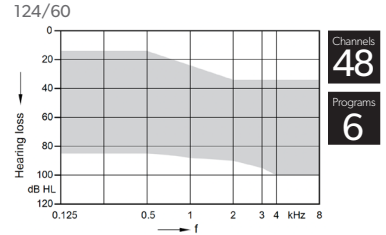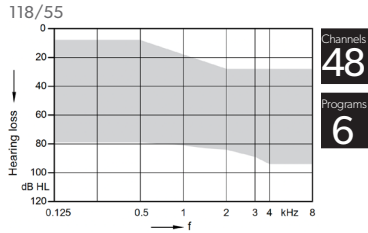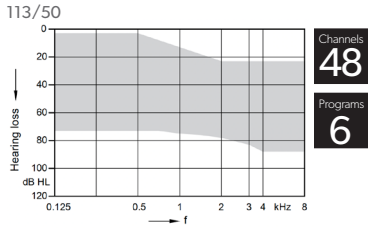

## IIC

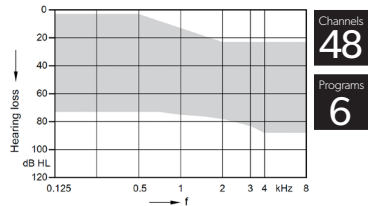

## Insigno Charge&Go AX Customs



Available in Full Shell, Half Shell, and Canal



**Insigno Charge&Go AX Charger (Standard)** – features contactless charging and 24 hours on a single charge



Black



Dark Brown



Brown



Mocha



Tan



Beige

### Key features & benefits



Rechargeable



iPhone & Android Connectivity



Immersive Soundscape



Tinnitus Relief



Augmented Speech Understanding



24/7 Signia Assistant



Automatic Situation Detection



Signia App



Auto EchoShield






HandsFree with CallControl for iOS



My WellBeing

Please note that red, blue and clear shells are also available.

## Inσιο Charge&Go AX Customs Fitting Ranges

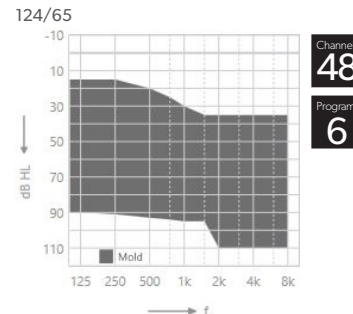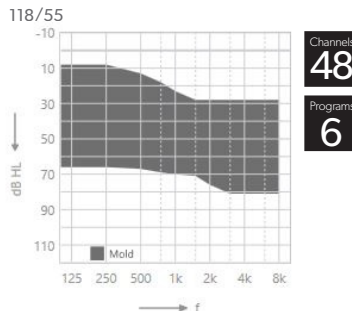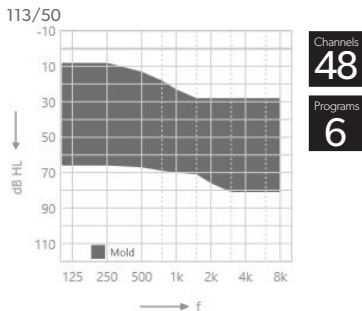

### Inσιο Indicators Available



Large R/L



Last Name



Large Red/Blue Dot\*



Canal Grips

### Inσιο Removal Options



Filament Loop



Retrieval Post



OtoRing

\*N/A on black and dark brown

# Insio Nx Custom



## Custom Faceplate Colors

Non-Bluetooth®  
Faceplate available

Description	ITE	ITC/HS	CIC	IIC
Beige	•	•		
Tan	•			
Mocha	•	•	•	
Brown	•	•	•	
Dark Brown	•	•		
Black				•

## Custom Shell Colors



Description	ITE	ITC/HS	CIC	IIC
Beige	•	•		
Tan	•		•	
Brown	•	•	•	
Dark Brown	•	•		
Clear	•	•	•	
Red	•	•	•	•
Blue	•	•	•	•

## Programmers



Stream audio directly into your Bluetooth® enabled Insio ITC/HS and ITE instruments with **Streamline™ TV** and **Streamline™ Mic**

## Insio Nx Custom Fitting Ranges

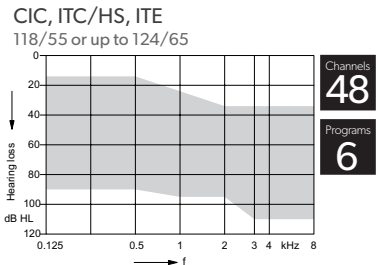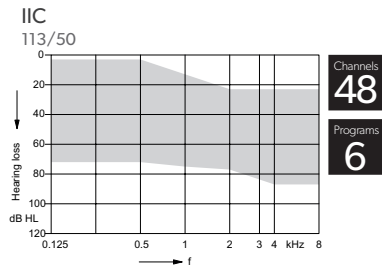







# Motion Charge&Go X

Complete BTE lineup with the world's first rechargeable Super Power

New EarWear 3.0 thin tubes, domes and sleeves are required for thin tube fittings



Standard

Power

Super Power



IP68



Signia Motion Charge&Go X



Signia Dry&Clean Charger for P & SP

Up to 61 hours of use on a single charge with SP, 30 hours for P and 24 hours for Standard (based on a full 4-hour charge)



Beige



Granite



Grey



Black



Silver Hero



Pearl White



Sandy Brown



Dark Champagne



Deep Brown



Rose Gold

## Key features & benefits



Lithium-Ion



iPhone & Android Connectivity



Integrated Telecoil



CROS



Tinnitus Relief



Dynamic Soundscape Processing



Acoustic Motion Sensors



Own Voice Processing



Signia App



Face Mask Mode



AI Digital Assistant


## Motion Charge&Go X Fitting Ranges

### Motion C&G X Fitting Range

60/133 dB SPL (Earhook)

60/125 dB SPL (Thin Tube 3.0)

63/126 dB SPL (ThinTube 3.0 P)

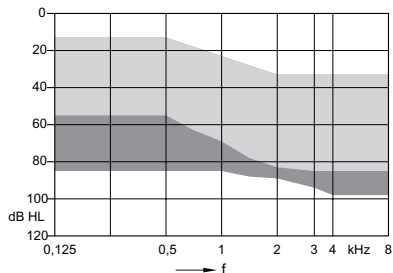


ThinTube 3.0 + Eartip 3.0 Open  
Earhook

### Motion C&G P X Fitting Range

77/135 dB SPL (Earhook)

66/130 dB SPL (ThinTube 3.0)


70/131 dB SPL (ThinTube 3.0 P)

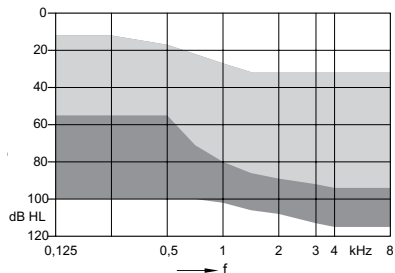

ThinTube 3.0 + Eartip 3.0 Open  
Earhook

### Motion C&G SP X Fitting Range

82/140 dB SPL (Earhook)

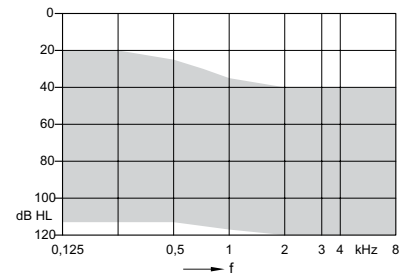

Earhook

## Motion 13 and 13 P Nx



IP68

Motion 13P Nx

Motion 13 Nx




Battery Type 13



Battery Type 13


### Fitting Ranges

70/139 – Motion 13 Nx

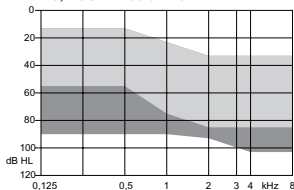

Channels **48**  
Programs **6**

81/141 – Motion 13P Nx

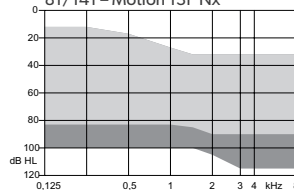

Channels **48**  
Programs **6**



Beige



Granite



Grey



Black



Silver



Pearl White



Golden Blonde



Sandy Brown



Dark Champagne




Dark Granite



Deep Brown

# BTE Earmolds

Bold text indicates default options



	BTE THIN TUBE 3.0 Earmold		BTE Conventional Earmold	
				
Styles	<b>Canal Hollow (Hard Material Only)</b> Canal Lock Hollow (Hard Material Only) Canal Canal Solid Helix Lock Solid	Full Shell Solid Half Shell Solid Semi-Skeleton Solid Skeleton Solid	Canal Solid Canal Lock Solid Helix Lock Solid Full Shell Solid	Half Shell Solid Semi-Skeleton Solid Skeleton Solid
Tubing	Thin Tube 3.0 (Size 0-4) Thin Tube 3.0 Power (Size 0-4)		<b>Standard</b> Thick Thick Dry Extra Thick Libby Horn	Standard Double Bend Clear Thick Double Bend Clear Thick Dry Double Bend Clear Extra Thick Double Bend Clear
Materials	<b>Hard</b>	Soft	<b>Hard</b>	Soft
Colors	<u>Gloss and Matte Finish Options</u> Clear, Rose, Beige, Brown, Red, Blue	<b>Clear Gloss is default</b>	<u>Gloss and Matte Finish Options</u> Clear, Rose, Beige, Brown, Red, Blue	<b>Clear Gloss is default</b>
Venting	None Open (Hollow Only) Pressure (1.0 MM) 1.2 MM 1.6 MM	SAV (2.0 MM) SAV (3.0 MM) Large as Possible Trench Vent (Hard Material Only) Semi-IROS	None <b>Pressure (1.0 MM)</b> 1.2 MM 1.6 MM SAV (2.0 MM)	SAV (3.0 MM) Large as Possible Trench (Hard Material Only) Semi-IROS
Ear Length	Short, <b>Medium</b> , Long (Hard Material Only)		Short, <b>Medium</b> , <b>Long</b>	
Indicators	<b>Red/Blue Dots</b> Red/Blue Letters	Patient Name	<b>Red/Blue Dots</b> Red/Blue Letters	Patient Name
Dexterity	Removal Line		Removal Line	
Wax Protection	<b>NanoCare 3.0</b>		None	

## Thin Tube 3.0

Compatible with Motion Charge&Go X and P BTE

Both ThinTube 3.0 options offer proven better retention, comfort, usability and physical fit rate

The new 1.4mm-width PowerTube also delivers a cosmetically appealing solution for more severe degrees of hearing loss.

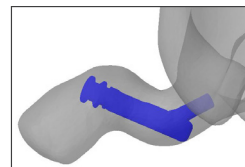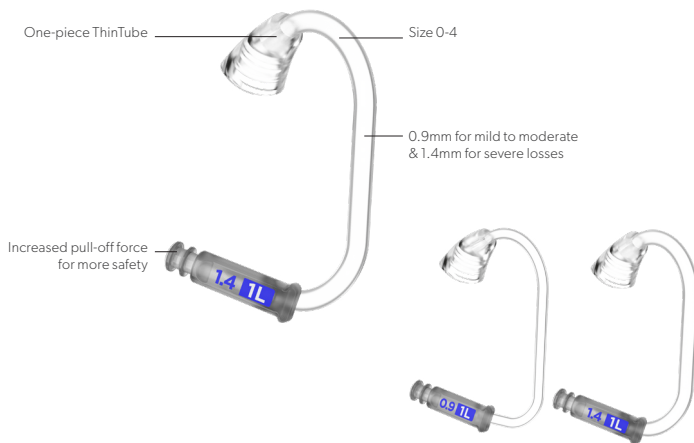

### Retention

Improved bending of tube + outlet locks the acoustical coupling reliably into the ear canal, prevents lateral migration and makes concha-locks redundant.

Binaural Thin tube 3.0 S0	10991609
Binaural Thin tube 3.0 S1	10991610
Binaural Thin tube 3.0 S2	10991611
Binaural Thin tube 3.0 S3	10991612
Binaural Thin tube 3.0 S4	10991613
Binaural Thin tube 3.0 S0 Power	10992492
Binaural Thin tube 3.0 S1 Power	10992493
Binaural Thin tube 3.0 S2 Power	10992494
Binaural Thin tube 3.0 S3 Power	10992495
Binaural Thin tube 3.0 S4 Power	10992496









## CROS



CROS Pure  
Charge&Go  
T IX/AX



CROS Pure  
312 AX



CROS  
Styletto IX



CROS  
Styletto AX



CROS Pure  
Charge&Go X  
compatible with  
Motion C&G X

- Li-ion rechargeability and auto on/auto off function when removed from charger for easy handling.
- 29 hours of battery life for CROS Charge&Go and 90 hours of battery life for CROS 312.
- World's first CROS/BiCROS technology to benefit from acoustic-motion sensors for the clearest speech understanding that adapts to the wearer's movements (with BiCROS).
- Own Voice Processing (OVP™) (with BiCROS).
- Directional pattern skewed toward the unaidable ear for improved speech understanding.
- CROS mixing-level adjustment available in the Signia App to balance spatial benefit vs. noise level.
- Compatible with Signia App, StreamLine™ TV, StreamLine™ Mic, and miniPocket™.

## CROS — RICs

DID YOU  
KNOW?

Signia leads the industry in CROS/BiCROS technology.

### A survey\* among Hearing Care Professionals found that:

- Signia CROS technology received an average of **4.5 out of 5 stars**, the majority stating it is the **BEST TECHNOLOGY** on the market.
- Most reported that they and their patients were **EXTREMELY SATISFIED** with Signia CROS.

### When asked for specific feedback, here's what your peers had to say:

- "It is **the only CROS product that I have had success with**. These patients still have some challenges, but overall they **accept the fitting at a much higher rate than any other CROS product** I have tried."
- "Patient is **hearing in situations they never heard** with their old [competitive CROS product]."
- Although we have used [competitor] in the past, my **previous [competitor] users who now have the Signia CROS are happier**.

\* Data on file.

## StreamLine Mic – Compatible with all Signia Bluetooth hearing aids

The new StreamLine™ Mic can be connected to virtually any Bluetooth®-enabled smartphone (such as iPhone® and Android™) and establishes various direct audio streaming options:

- Phone call streaming (both ears)
- Integrated Microphone
- Hands-free phone calls (both ears)
- Accept/decline calls
- Remote microphone functionality (both ears)
- Hear in challenging situations with remote microphone mode
- Audio, video and music streaming (both ears, stereo)
- Remote Control functions
- Volume Increase
- Volume Decrease
- Program Change

“Made for iPhone” means that an electronic accessory has been designed to connect specifically to iPhone, and has been certified by the developer to meet Apple performance standards. Apple is not responsible for the operation of this device or its compliance with safety and regulatory standards. Please note that the use of this accessory with iPhone may affect wireless performance. iPhone is a trademark of Apple Inc., registered in the U.S. and other countries.

The Android name, the logo, and other trademarks are property of Google LLC.



Watch handling videos online  
[veterans.signiausa.com/howto](https://veterans.signiausa.com/howto)

### HOW TO ORDER IN ROES

ROES Main Menu – Patient Orders.

Select an Action: Wireless Devices

Remotes

Wireless Transmitter

Select Device Type: SIVANTOS

Select Make: StreamLine Mic

## StreamLine Mic – Overview

### Multi-function button (MFB)

- Switch on (press 3 sec) and off (press 6 sec)
- Toggle between two modes:
  - Bluetooth® Phone Mode
  - Remote Microphone Mode


### Wearer controls/status indicators

- VC / Program up
- VC / Program down
- Power indicator
- Status indicator






### StreamLine Mic offers intuitive onboard control and status elements
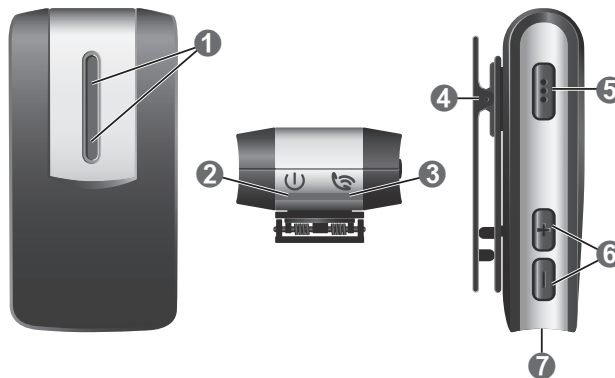
### Accessories included in the shipment

- 2 lanyards (regular and long)
- USB wall plug + USB connector
- Pouch
- User Guide



### StreamLine Mic controls

- 1 Microphones
- 2 Power indicator 
- 3 Status indicator 
- 4 Lapel clip
- 5 Multi-function button (volume, program, mute) 
- 6 Volume buttons  and 
- 7 USB port for charging

## StreamLine Mic – Pairing Quick Start Guide

### Pairing with Hearing Aids

- Open/close the battery doors of the hearing aids and/or place in/out of the charger to place them in pairing mode.
- Turn on StreamLine Mic by holding down the MFB until the power indicator is green. This takes approximately 3 seconds. The status indicator flashes yellow.
- Place both hearing aids close to StreamLine Mic (maximum distance 20 cm or 8 inches). The pairing process starts automatically. Pairing is finished when the status indicator is solid yellow.
- A solid green power indicator and solid yellow status indicator verify connection to hearing aids is successful with no mobile device paired.
- Your hearing aids need to be paired once with StreamLine Mic.

### Pairing with Smartphone

- Confirm that Bluetooth® of the mobile device is on.
- Turn on StreamLine Mic by holding down the MFB until the power indicator is green. This takes approximately 3 seconds.
- Hold down the the MFB and “+” button on the StreamLine Mic at the same time until status indicator flashes blue. This takes approximately 6 seconds.
- The wearer has 3 minutes to complete the pairing.
- Select “StreamLine Mic” from the list of detected devices on the Bluetooth® menu.
- If requested, enter the PIN “0000”.
- A solid green power indicator and solid blue status indicator verify connection to smartphone is successful.



### StreamLine Mic – Battery Runtime

- Charging status is clearly displayed via power indicator on the StreamLine Mic
- Full charge takes approx. 2 hours
- Up to 6 hours battery runtime after full charge with full functionality (phone calls, audio streaming, remote mic)
- Can be charged via USB or usual plug-in (all necessary cables are included in the delivery)
- Can be used while charging

## StreamLine TV – Compatible with all Signia Bluetooth® hearing aids

The perfect companion for TV fans.

### StreamLine TV.

Thanks to StreamLine™ TV audio from television or any device with audio output can be fed directly into Signia hearing aids. Wearers can even adjust the volume individually via the Signia app.

All wearers of IX, AX, X and Nx hearing aids have direct binaural streaming for TV via the StreamLine TV.

This includes use without the Signia app for non-smartphone users.



#### HOW TO ORDER IN ROES

ROES Main Menu – Patient Orders.

Select an Action: Wireless Devices

Remotes

Wireless Transmitter

Select Device Type: SIVANTOS

Select Make: StreamLine TV



## Pairing for StreamLine TV

If ordered together, Signia will pair all Bluetooth® compatible hearing aids to the Streamline TV.


The hearing aids need to be paired once with the transmitter. After that, the pairing is saved in the hearing aids and remembered each time they are turned on.

- Ensure that the transmitter is connected to the power supply.
- Turn the hearing aids off and on again: open and close the battery compartments. If your hearing aids do not have a battery compartment, turn them off and on as described in the hearing aids' user guide.
- Place both hearing aids close to, or on top of, the transmitter. The pairing process starts automatically.
- Observe the indicators:



### Description

- ■ If pairing is completed, both indicators are blue for about 30 seconds.

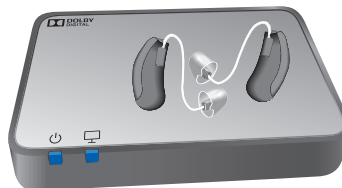


### Multiple users:

- The transmitter can be paired to multiple hearing aids. All hearing aids will be able to hear the same audio signal.
- To clear the pairings of one current user, first reset the transmitter to factory settings. This will clear the pairings of all users. After that, pair the hearing aids of the other users again. Refer to the section "Resetting to factory settings".

## Connecting StreamLine TV

### TOSLINK connection

- Before connecting the cable, remove the protection caps from the TOSLINK cable plugs.
- Plug the cable into the TOSLINK **input** of the transmitter and into the TOSLINK (optical audio) **output** of electronic device. To prevent the TOSLINK cable from breaking, do not forcefully bend it.

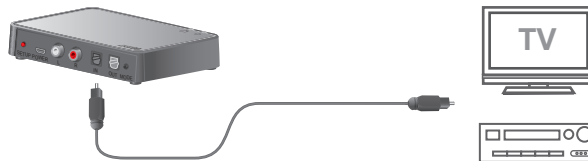


### Further information

- Refer to the electronic device's manual: set the output format of the electronic device to either PCM (pulse-code modulation) or Dolby® Digital.\*

### RCA connection

- Plug the cable into the RCA connections of the transmitter.
- Plug the other side of the cable into the RCA output of the electronic device or via adapter into the SCART or via mini jack output of the electronic device.

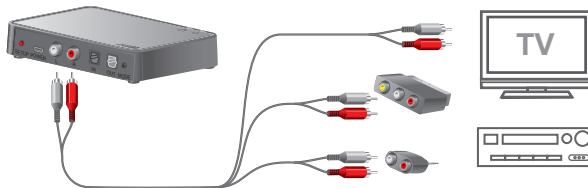

For more information and details on how to connect to your home theater amplifier or a sound bar watch the handling videos online at [veterans.signiausa.com/howto](https://veterans.signiausa.com/howto)

\*Dolby and the double-D symbol are registered trademarks of Dolby Laboratories.

## miniPocket — Remote control compatible with all Signia hearing aids


Small enough to hang on a keyring, the miniPocket™ is ideal for wearers who want quick, discreet control of their hearing aids without using a smartphone.

- Battery type CR2450
- Pairing remains in place when battery dies and is replaced

---

### Components

- 1 Status LED
  - 2 Volume buttons
  - 3 Keyring
  - 4 Key lock switch
  - 5 Program change button
- 

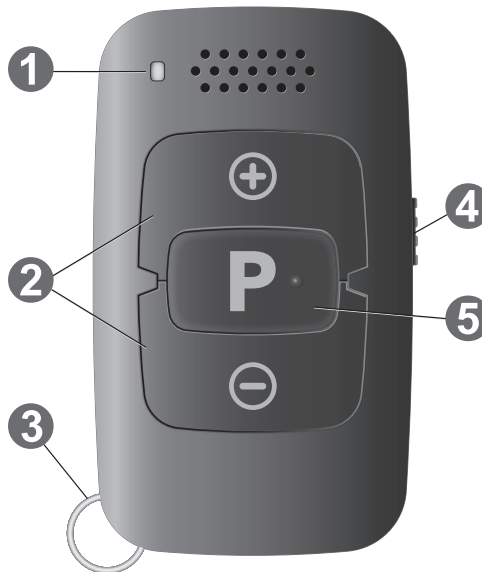

## Pairing hearing aids with miniPocket

1. Move key lock switch to 'lock' position



2. Turn on hearing aids and place in ears
3. Hold down **P** button on the remote control while switching the key lock to the unlock position
4. Release **P** button (LED on)
5. Press and hold **P** button until the pairing confirmation tone is heard
6. Move Key Lock switch to Lock position to exit pairing mode



Watch video online

[veterans.signiausa.com/howto](https://veterans.signiausa.com/howto)





### HOW TO ORDER IN ROES

When ordering a hearing aid in ROES, look for the "Select a remote" dropdown. Select "miniPocket"

Chargers



## Charger Compatibility

Chargers		Integrated Xperience (IX) and Augmented Xperience (AX)			
		Pure Charge&Go Portable Charger*	Pure Charge&Go Charger	Pure Charge&Go Dry&Clean Charger	Insiio Charge&Go AX Charger
Type of charger		Galvanic	Galvanic	Galvanic	MR Inductive
Part number		10993118	10996345	10991650	21024145
					
Cord connection		USB-C PN 21023081	Attached	USB-C PN 21023081	USB-C PN 21023081
Wear time without streaming		39h	39h	39h	24h
Wear time with streaming		36h	36h	36h	20h
Charge time		4h	4h	4h	4h
Quick charge time		30m for 9h	30m for 9h	30m for 9h	30m for 4h
Portable charger extra charges		3	N/A	N/A	N/A
Hearing Aids					
RIC	Pure C&G IX/AX	*	•	•	
	CROS Pure C&G IX/AX	*	•	•	
	Styletto AX				
	CROS Styletto AX				
Custom	Insiio C&G AX				**

\* Compatible \* Default (Recommended)

Please note – While the X hearing aids will work with the Inductive Charger I, the charger was not custom made for their shape. Please verify that the lights turn on when placing the hearing aid in the charger so that you can ensure proper placement.

## Charger Compatibility

Chargers					
	<b>Styletto IX Charger</b>	<b>Styletto AX Charger*</b>	<b>Styletto AX Dry&amp;Clean Charger</b>	<b>Charging Pad<sup>1</sup></b>	<b>Clinic Multicharger</b>
Type of charger	Galvanic	Galvanic	Galvanic	-	Galvanic
Part number	10997503	21003216	21007593	PN 10245297	PN 21026076



Cord connection	USB-C PN 21023081	USB-C PN 21023081	USB-C PN 21023081		Attached
Wear time without streaming	24h	20h	20h		39h
Wear time with streaming	20h	17h	17h		36h
Charge time	4h	4h	4h		4h
Quick charge time	30m for 4h	30m for 5.5h	30m for 5.5h		30m for 9h
Portable charger extra charges	5	4	N/A		N/A

Hearing Aids						
RIC	Pure C&G IX/AX					•
	CROS Pure C&G IX/AX					•
	Styletto IX	•		•		
	Styletto AX		*	•	•	•
	CROS Styletto AX		*	•	•	•

\* Compatible    <sup>1</sup> Default (Recommended)

Please note – While the X hearing aids will work with the Inductive Charger I, the charger was not custom made for their shape. Please verify that the lights turn on when placing the hearing aid in the charger so that you can ensure proper placement.



## Charger Compatibility

Chargers	Xperience (X)								
	Inductive Charger I	Inductive Charger II	Styletto X Charger	Dry&Clean Charger P X	Motion C&G Charger P X	Dry&Clean Charger SP X	Motion C&G Charger SP X	Charging Pad <sup>1</sup>	Clinic Multicharger
Type of charger	Inductive	Inductive	Galvanic	Galvanic	Galvanic	Galvanic	Galvanic	-	Galvanic
Part number	10948199	10966942	10972015	10991638	10996335	10991644	10996340	PN 10245297	PN 21026076
Cord connection	micro USB	micro USB	USB-C	USB-C	Attached	USB-C	Attached		Attached
Wear time without streaming		24h	19h	30h	30h	61h	61h		39h
Wear time with streaming		21h	16h	27h	27h	57h	57h		36h
Charge time		4h	4h	4h	4h	4h	4h		4h
Quick charge time		30m for 6h	30m for 5h	30m for 7h	30m for 7h	30m for 4h	30m for 4h		30m for 9h
Portable charger extra charges		N/A	3	N/A	N/A	N/A	N/A		N/A
Hearing Aids									
RIC	Pure C&G X	•	••						
	CROS Pure C&G X	•	••						
	Styletto X			••				•	
	Pure C&G Nx	•	••						
BTE	Motion C&G X	•	*						
	Motion C&G P X			••	•				•
	Motion C&G SP X					••	•		•

• Compatible • Default (Recommended)

Please note – While the X hearing aids will work with the Inductive Charger I, the charger was not custom made for their shape. Please verify that the lights turn on when placing the hearing aid in the charger so that you can ensure proper placement.

## Rechargeability



### Pure Charge&Go T IX/AX Standard Charger

- Compatible with Pure Charge&Go T IX, Pure Charge&Go T AX and CROS Charge&Go AX
- Convenient, table-top design



### Pure Charge&Go T IX/AX Portable Charger

- Compatible with Pure Charge&Go T IX, Pure Charge&Go T AX and CROS Charge&Go AX
- Portable, closed-lid charger
- Holds three charges



### Insio Charge&Go AX Charger

- Compatible with Insio Charge&Go AX featuring contactless charging



### Inductive Charger II

- Compatible with Pure Charge&Go X and Motion Charge&Go X
- Electronically dehumidifies the hearing aids
- Room for earmolds



### Styletto IX/AX Charger

- Compatible with Styletto IX/AX
- Portable, closed-lid charger equipped with Qi Wireless Technology — enables wireless charging on any standard charging pad, including those built into cars and cafe counters\*



### Dry&Clean Charger

- Dry&Clean Chargers are available for Motion C&G P, Motion C&G SP, and Pure C&G T AX. This charger will simultaneously charge, electronically dehumidify, and sanitize the hearing aids and custom molds with the integrated UV-C light.
- Compatible with Motion X, Pure IX, Pure AX, and Styletto AX

\*Wireless charging pad not included with purchase

## Portable IX/AX Charger

For Pure Charge&Go T IX, Pure Charge&Go T AX, and CROS Charge&Go IX/AX



### NOTICE

Please make sure that the power plug is easily accessible to remove it from power supply if necessary.

- Connect the charger to the power supply. All LEDs for each charging slot turn on successively indicate the charger has power and is ready for use.
- Put the hearing aids in the charging slots as shown in the picture. You can charge both hearing aids at the same time. The hearing aids are automatically turned off and start charging. The charging status is indicated, see table "Charging status".
- When you want to use your hearing aids, take them out of the charger. They will automatically turn on.



## Charger Lights

---

### LED when charger is charging:

---

**Pulsing green:** internal battery is being charged.

---

**Solid green:** charger is fully charged.

---

**Blinking red:** charger error. Please refer to section "Troubleshooting".

---

---

### LED when you are mobile:

---

**Solid green:** the internal battery has enough capacity to fully charge a pair of hearing aids.

---

**Solid orange:** the internal battery does not have enough capacity to fully charge a pair of hearing aids.

---

**Blinking orange:** the internal charger battery is empty. Charging of hearing aids is not possible.

---

## Inσιο Charge&Go Charger

For Inσιο Charge&Go AX












### NOTICE

Please make sure that the power plug is easily accessible to remove it from power supply if necessary.

- Connect the charger to the power supply. All LEDs for each charging slot turn on successively indicate the charger has power and is ready for use.
- Put the hearing aids in the charging slots as shown in the picture. You can charge both hearing aids at the same time. The hearing aids are automatically turned off and start charging. The charging status is indicated, see table "Charging status".
- When you want to use your hearing aids, take them out of the charger. They will automatically turn on.



## Charging status

LED	Status description
	<b>No LED:</b> Charger is not connected to power supply.
	Single <b>green</b> LED turns on successively: Charger is connected to a power supply and ready to charge. No hearing aid is inserted in the respective charging slot or the hearing aid is not inserted correctly.
	Flashing <b>orange</b> LED: Charging status of the hearing aid is read out before charging starts.
	Single <b>orange</b> LED: Charging status 0 to 33 %
	Turn on single <b>green</b> LED: Charging status 34 to 66 %
	Turn on two <b>green</b> LEDs: Charging status 67 to 99 %
	Three <b>green</b> LEDs: Hearing aid is fully charged.
	Single <b>orange</b> LED flashes successively: Cooling phase Charging is resumed after 5 minutes.
	Flashing <b>red</b> LEDs: Error Read more in section "Troubleshooting".

## Inductive Charger II

For Pure Charge&Go X, CROS Charge&Go X, Motion Charge&Go X BTE



### NOTICE

Please make sure that the power plug is easily accessible to remove it from power supply if necessary.

- Connect the charger to the power supply. All LEDs for each charging slot turn on successively indicate the charger has power and is ready for use.




- Put the hearing aids in the charging slots as shown in the picture. You can charge both hearing aids at the same time. The hearing aids are automatically turned off and start charging. The charging status is indicated, see table "Charging status".
- When you want to use your hearing aids, take them out of the charger. They will automatically turn on.

## Charging status

LED	Status description
	<b>No LED:</b> Charger is not connected to power supply.
	Single <b>green</b> LED turns on successively: Charger is connected to a power supply and ready to charge. No hearing aid is inserted in the respective charging slot or the hearing aid is not inserted correctly.
	Flashing <b>orange</b> LED: Charging status of the hearing aid is read out before charging starts.
	Single <b>orange</b> LED: Charging status 0 to 33 %
	Turn on single <b>green</b> LED: Charging status 34 to 66 %
	Turn on two <b>green</b> LEDs: Charging status 67 to 99 %
	Three <b>green</b> LEDs: Hearing aid is fully charged.
	Single <b>orange</b> LED flashes successively: Cooling phase Charging is resumed after 5 minutes.
	Flashing <b>red</b> LEDs: Error Read more in section "Troubleshooting".

## Styletto IX/AX Charger

- Place the right hearing aid in the charging slot marked with an "R", and the left hearing aid in the charging slot marked with an "L".
- Stow the ear pieces nicely.
- Close the charger's lid for battery power saving reasons.
- The hearing aids are automatically turned off and start charging. The charging status is shown for 7 seconds.

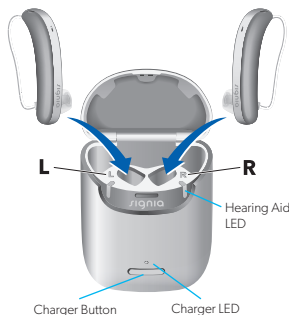

If the status LEDs do not light up when you insert the hearing aids, check the following:

- The hearing aids are not correctly placed. Push them slightly into their charging slots or close the charger's lid.
- The charger is turned off. Press and hold the button for 5 seconds to turn on the charger.
- The charger's internal battery is too low for mobile charging. Connect the charger to an electrical outlet.



**NOTICE:** If you put a hearing aid in the wrong charging slot, the receiver cable may be damaged when closing the charger's lid.

### Charging status of hearing aids

The LED next to the hearing aid indicates the charging status of the hearing aid. Note that the status is only displayed immediately after putting the hearing aids in the charger, when you briefly press the button or when you open the lid. When the charger is off, no status is displayed at all.

#### Hearing Aid LED Status description

**Pulsing green:** Hearing aid is being charged.

**Green:** Hearing aid is fully charged. Internal battery has enough capacity to fully charge two hearing aids.

### Charger LED Status description

**LED is off:** No hearing aid is in the charger, hearing aid is not placed correctly, or charger is off.

**Pulsing orange:** Cooling phase. Charging is resumed when the device has cooled down. Please refer to section "Troubleshooting".

**Blinking red:** Hearing aid error. Please refer to section "Troubleshooting".

**Green:** Internal battery has enough capacity to fully charge two hearing aids.

## Dry&Clean Charger

For Pure Charge&Go T IX, Pure Charge&Go AX, Styletto IX/AX, Motion P X and Motion SP X



### Cleaning and Drying

- Drying and cleaning focuses on tubes, ear tips and molds.
- Put the hearing aids into the charging slots and close the lid.
- UVC cleaning starts automatically after closing the lid and lasts 15 minutes.
- Ongoing UVC cleaning process is indicated by LEDs (see table "Charging and cleaning status").
- Drying starts automatically after closing the lid and lasts 4 hours.


### Charging

A full charging cycle takes approximately 3 to 4 hours.





After 30 minutes, it is possible to use your hearing aids for up to 7 hours.

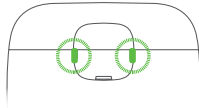


Charging times may vary due to the use of the hearing aids and the age of the power cell.


## Dry&Clean Charger






### CHARGE INDICATORS

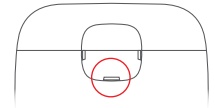

<p>Solid yellow</p> 	The charger is detecting the hearing aid.
<p>Flashing green</p> 	The hearing aid is charging, low battery, less than 40%.
<p>Glowing green</p> 	The hearing aid is charging, battery more than 40%.
<p>Solid green</p> 	The hearing aid is fully charged, battery 100% charged.

<p>One indicator flashing red</p> 	An error occurred during charging, and the charging has stopped.
<p>Both indicators flashing red</p> 	An error occurred during charging, and the charging has stopped.
<p>Flashing yellow</p> 	The temperature of the hearing aid has become too high during charging, and charging is paused.



### DRY AND CLEAN INDICATORS

<p>Solid yellow</p> 	The maintenance cycle is paused, but the charging cycle continues.
<p>Flashing white</p> 	The maintenance cycle is in progress (drying/UV-C).
<p>Solid white</p> 	The maintenance cycle is complete.

<p>Flashing red</p> 	An error occurred during the maintenance, the cycle has stopped.
<p>Flashing yellow</p> 	The maintenance cycle is paused due to overheating (fan still running).



## How to order a Charger in ROES

Choose from Other Commodities, and follow these ordering steps:

ROES Main Menu – Patient Orders

**Select an Action:** Other Commodities or Accessories

**In Commodity Group pull-down menu, choose:** SIVANTOS

**In Vendor pull-down menu, choose:**

In Search box please enter one of the following words:

ROES Name	Description	Part Number
Insio Charge&Go AX Charger	Charger for Insio Charge&Go 7 AX	21025145
Pure Charge&Go IX/AX Portable Charger	Compatible with Pure Charge&Go T IX/AX, CROS Pure IX/AX	10993118
Pure Standard Charger	Charger for Pure Charge&Go T IX/AX, CROS Charge&Go IX/AX	10996345
Pure Dry&Clean Charger	Charger for Pure Charge&Go T IX/AX, CROS Charge&Go IX/AX	10991650
Styletto Charger	Charger for Styletto IX/AX, CROS Styletto IX/AX	10972015
Inductive Charger II	Charger for Pure Charge&Go X, CROS Pure X, and Motion Charge&Go X	10966942
Dry & Clean Charger Motion X BTE P	Dry & Clean Charger for Motion Charge&Go P X BTE	10991638
Dry & Clean Charger Motion X BTE SP	Dry & Clean Charger For Motion Charge&Go SP X BTE	10991644
Smart Lithium-ion Dryer	NX and older charging options	10960141

**NOTE** If you want it delivered to your station, then change the “Deliver To” shipping address, otherwise, it will ship to the patient by default. Please note: a charging unit comes shipped with every Li-ion hearing aid order.

## Portable Power Supply

- The Portable 2600mAh Mobile Power Bank is the perfect power solution when on the road or stuck somewhere without convenient electricity.
- It fits comfortably in a pocket, purse, backpack, or laptop bag to keep Pure Charge&Go or Motion Charge&Go hearing aids charged when needed.
- Capacity of 2600mAh will provide three full charges.



Portable Power Supply (10948100)








## Signia Nomenclature

The names of our hearing aids are clues to what's under the hood....

Make	Style	Model	Feature
Pure	RIC	C&G IX T	Rechargeable Integrated Xperience Telecoil
Stylecto	Slim RIC	IX AX	Rechargeable Integrated Xperience Rechargeable Augmented Xperience
Insio	Custom	IX IIC/CIC C&G AX	Integrated Xperience Rechargeable Augmented Xperience
Motion	BTE, P BTE, SP BTE	C&G X	Rechargeable Xperience
Pure	RIC	312 AX	Battery AX Platform
Insio	Custom	NX	Battery NX platform
CROS	RIC Custom	C&G, 312, Stylecto C&G IX	IX transmits to IX AX transmits to AX X transmits to X C&G transmits to C&G 312 transmits to 312


## Where to Find Serial Numbers

C&G RIC	Under the Rocker Switch
RIC 312	Under the Rocker Switch
Customs	Printed on shell
Styletto AX	Underside of aid
Streamline™ Mic	Under swivel clip
Streamline™ TV	Under SLTV
MiniPocket™	Sticker in battery compartment

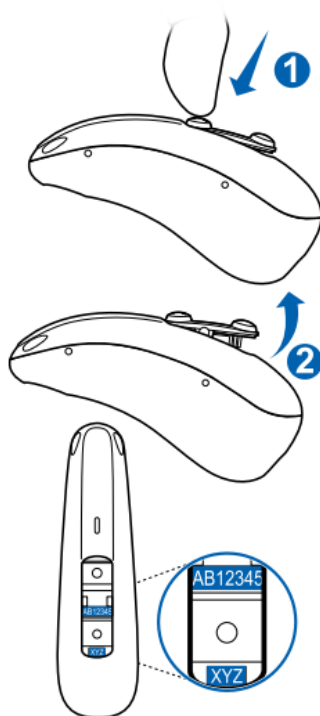





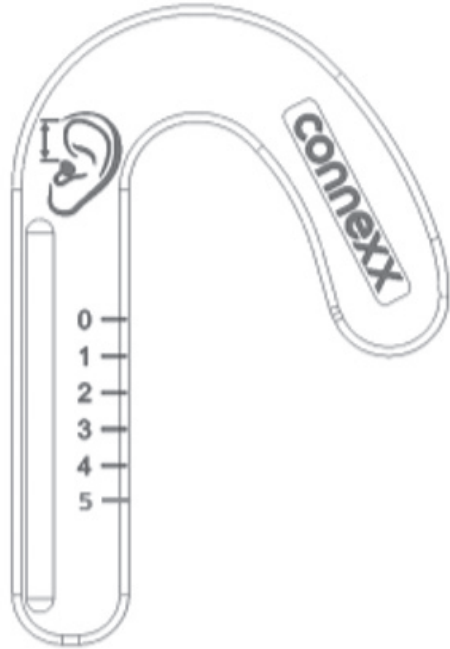

## How to Measure for Signia 3.0 Receivers

1. Use Measuring Gauge for EarWear 3.0 (white)
2. Place the measuring gauge on the ear of the wearer.
3. Read the number that comes closest to the upper portion of the ear canal opening.



## Signia Earwear 3.0 Parts

Item	Description	Brand		
Measurement Gauge 3.0		Ensure perfect match for new miniReceiver/Thin Tube 3.0 in 5 lengths (0-4)	Connexx	
Wax Guard 3.0 NanoCare – miniReceiver 3.0			Matchbook design containing 8 wax guards	Connexx, Signia
CeruSTOP 3.0 XL – for HP Earmolds		10 pieces per bag	Connexx	
Removal Tool 3.0		Tool to remove miniReceiver/Thin Tube from Earmold 3.0 Short/Long		
Handling Tool 3.0 for custom eartip		Tool to press out the receiver from Earmold 3.0 Short		



## Signia Wax Guard Guide



**miniR Wax Guard**  
Part #1096257  
Used for: 2.0 Receivers




**3.0 NanoCare**  
Part #10993607  
Used for: 3.0 Receivers and 3.0 Earmolds



**QuickGuard**  
Part #10962504  
Used for: Insignia Customs and 2.0 Earmolds




**3.0 CeruSTOP XL**  
Part #10999334  
Used for: 3.0 HP Receivers

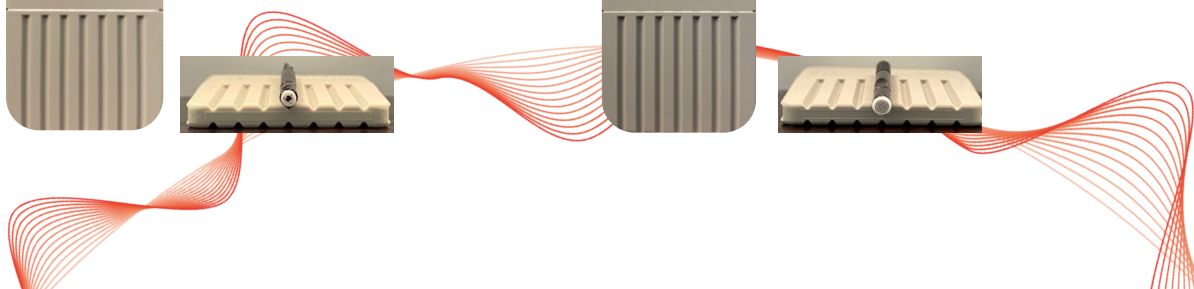


## Earmold Ordering Tips

- Order 3.0 Earmolds for all IX and AX receivers and Motion X Thin Tubes
- 3.0 Hollow Earmolds are open earmolds (same as Foil Earmolds for the 2.0 receivers)
- All HP receivers need an impression as they have to have earmold
- Our CROS has a real receiver and any earmold will work. We suggest Hollow for 3.0 earmolds or Foil for the 2.0 earmolds

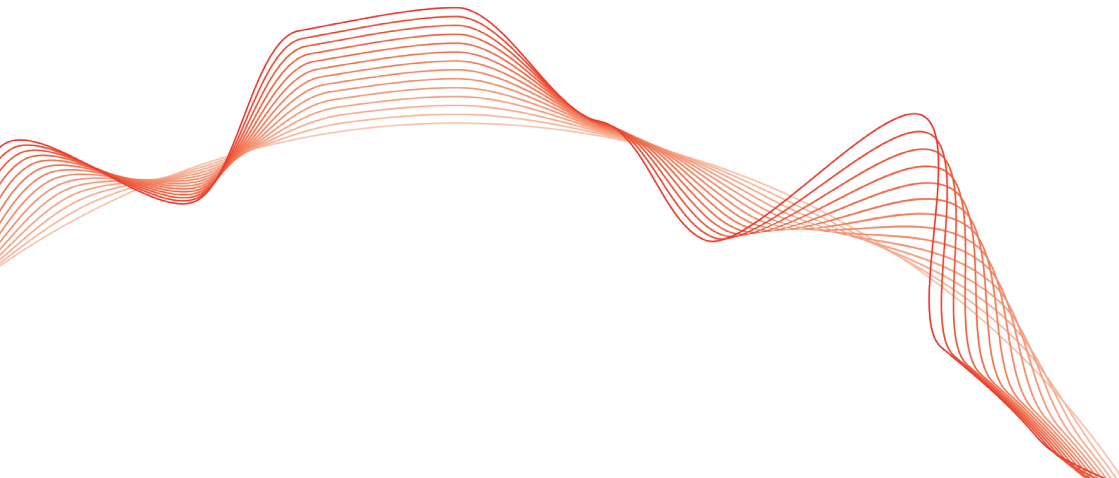

## Accessories Tips

- The miniPocket takes a CR2450 battery – this is labeled as an “L2450” battery in ROES for re-ordering
- Current accessories do not have to be connected/paired to the aids in the software
- The medium press of the hearing aid push button/rocker switch is defaulted to TV Stream for easy access

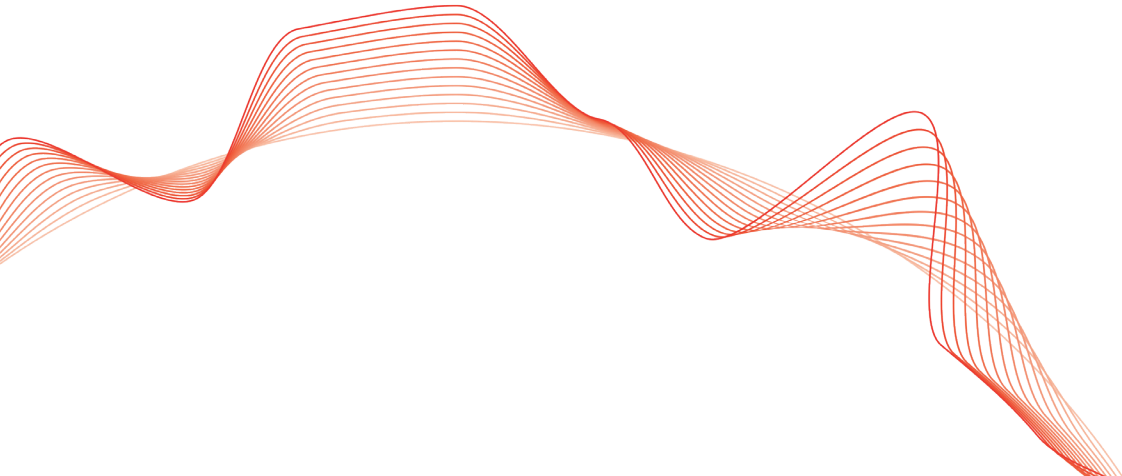

# Signia Fitting Guide

**NOTE: All ordered hearing aids have already been First Fit.**

1. Select your patient from NOAH patient database
2. Launch Signia Connex fitting software
3. Select "Detect Hearing aids"



Audiogram




Detect Hearing Instruments




Hearing Loss Simulator

4. Once hearing aids are detected, read out connected devices

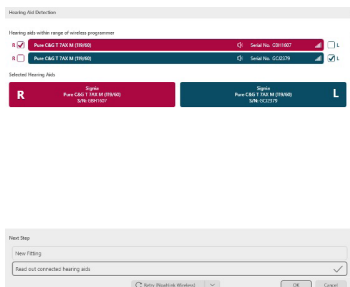


Please note:

- a. Since aids have been first fit, if you are reading out connected instruments (recommended) Connex will open on the Data Log screen.
- b. Visually confirm green check mark beside First fit.

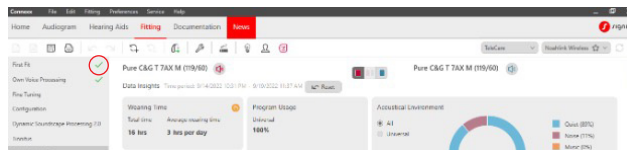

\*Note: Bone conduction is considered for fitting. If Audiogram in Noah does not show bone conduction, then you may get an audiogram error message.

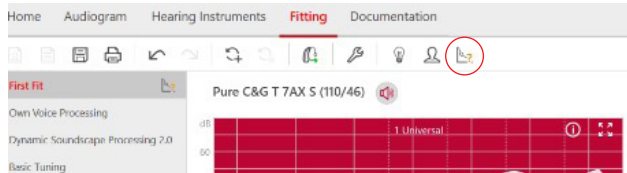



- c. If any acoustic parameters are adjusted on the First Fit tab, we recommend performing a new first fit.

# Signia Fitting Guide


(continued)

## 5. Run OVP (Own Voice Processing)\*

For optimal results

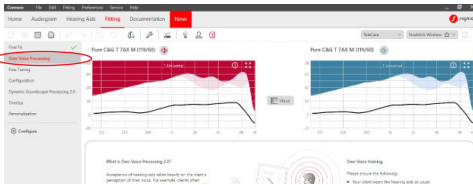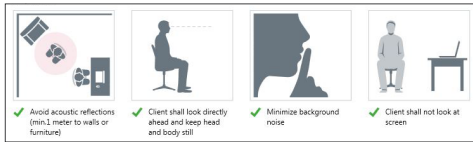

Complete the OVP Training. Tell your patient that they will hear a beep and that they should count loudly and continuous from 21 upwards. Example: 21, 22, 23, and so on.

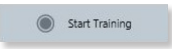

OVP is automatically enabled. Select "OK". Training was successful

\*OVP is not available with customs.

## 6. Recommended Best Practices for adjusting gain to Real Ear Targets

a. While aids are coupled adjust overall gain for aids in Fine Tuning > Sound & Loudness, Frequency Shaping tab under gain curve screen.



b. If further fine tuning is needed uncouple Right and Left said go into Gain tab under gain curve and adjust soft moderate and loud gain curves independently.




# Signia Fitting Guide

(continued)

## 7. Select Configuration Configuration

- **Hearing aid** – Configure User Controls, Volume Control, Therapy Signal, SoundBalance.
- **System Sounds** – Configure Alert Tones.

8. Save Session and Select “X” at top right of screen or select drop down item under File.
  - Disconnect hearing aids.

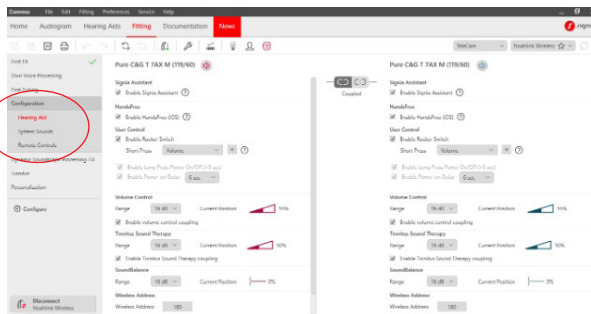

# Tinnitus Fitting Guide

Notch Therapy is only available in Signia hearing instruments. Notch Therapy can be customized to treat the patient's tonal tinnitus in a few simple steps with Connexx® fitting software. Find the patient's pitch match by using Guided Matching, Manual Matching, or choosing a known pitch from Direct Entry. Perform the frequency check to verify the correct octave of the pitch match and finally apply the notch to desired programs. It's that quick and easy.

The screenshot shows the 'Fitting' tab in the Connexx software. The main display area shows a graph titled 'Pure C&G T 7AX M (119/60)' with a notch at 3333 Hz. The graph plots dB HL (0 to 120) against frequency (125 to 8k Hz). Below the graph, the 'Guided Matching' section shows a 'Completed' status: 'Frequency 3333 Hz has been recorded as tonal tinnitus.' To the right, the 'Frequency Check' section asks 'Which tone sounds most similar to your perceived tinnitus?' and lists three tones: Tone A (35 / 35 dB HL, 1667 Hz), Tone B (38 / 38 dB HL, 3333 Hz), and Tone C (36 / 36 dB HL, 6667 Hz). The 'Activate Notch' section is also visible.

## CROS Quick Fit Guide

Under Fine Tuning, select Sound Settings to adjust the sensitivity of the CROS microphone in relation to the hearing aid mic.

If the “Enable CROS microphone” is unchecked, the CROS microphone will be disabled.

### Input Configuration

Mode: CROS

Enable CROS microphone

CROS level







## CROS/BICROS Listening Check

### Outside of Connexx

1. Turn on Both Hearing aid and CROS/BICROS
2. Hearing Aids need to be “Head Distance Apart”
3. Rub CROS/BICROS Microphones

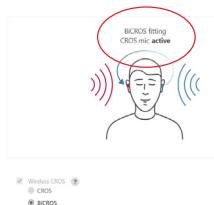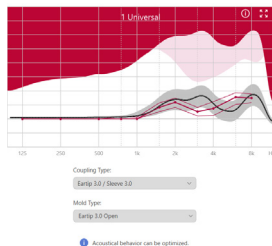

### Inside Connexx Software

1. Detect/Connect to both Hearing Aid and CROS/BICROS
2. First Fit Settings will show CROS/BICROS Mic Active
3. Fine Tuning>Sound Settings > CROS Mic Slider Enabled

### Input-Configuration

Mode: CROS

Enable CROS microphone

 CROS Microphone is not available in

CROS level

-12 dB

Note: Start up delay defaults to 18 seconds.

## How to Pair

### Pairing with an iPhone

1. Open the “Settings” app on the iPhone, OR from the Signia app select “Bluetooth®” connection option.
2. In the settings menu on your phone:
  - a. Select “Accessibility”
  - b. Select “MFi Hearing Devices” or “Hearing Devices”
3. Ensure Bluetooth® is enabled and pair the hearing aids when prompted; the names of your hearing aids will be in blue when successful.
4. Choose desired connection method.
5. Open the Signia app and tap “Let’s Go” to pair your hearing aids to the app.

### Pairing with an Android

1. When pairing your hearing aids to an Android device, go directly into the Signia app (the phone’s native Bluetooth® menu does not allow correct pairing).
2. Open Signia app and select “Bluetooth®” connection option.
3. Verify phone is on the list, tap “Next”.
4. Restart the hearing aids and place them next to the phone. Tap “OK”.
5. Tap the names of your hearing aids once they appear.
6. Tap “Let’s Go” to pair your hearing aids to the app.



Visit [signiausa.com/signia-app](https://signiausa.com/signia-app) for interactive user guide and additional information:







## ROES codes

To place an order for any of the following in ROES first choose **Commodities**, then **Parts and Accessories**, then search by item number.

VA ROES ITEM	DESCRIPTION	ITEM #	QTY.	PRICE	COMPATIBLE PRODUCTS
<b>RECEIVERS</b>					
MINI RECEIVER 3.0 KIT	MINI RECEIVER KIT 3.0 Complete	21018572	1	\$728.00	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 45 dB SIZE 0 RIGHT EAR	MINIRECEIVER 3.0 45 dB SIZE 0 RIGHT EAR	10992548	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 45 dB SIZE 0 LEFT EAR	MINIRECEIVER 3.0 45 dB SIZE 0 LEFT EAR	10992549	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 45 dB SIZE 1 RIGHT EAR	MINIRECEIVER 3.0 45 dB SIZE 1 RIGHT EAR	10992550	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 45 dB SIZE 1 LEFT EAR	MINIRECEIVER 3.0 45 dB SIZE 1 LEFT EAR	10992552	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 45 dB SIZE 2 RIGHT EAR	MINIRECEIVER 3.0 45 dB SIZE 2 RIGHT EAR	10992553	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 45 dB SIZE 2 LEFT EAR	MINIRECEIVER 3.0 45 dB SIZE 2 LEFT EAR	10992554	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 45 dB SIZE 3 RIGHT EAR	MINIRECEIVER 3.0 45 dB SIZE 3 RIGHT EAR	10992555	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 45 dB SIZE 3 LEFT EAR	MINIRECEIVER 3.0 45 dB SIZE 3 LEFT EAR	10992556	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 45 dB SIZE 4 RIGHT EAR	MINIRECEIVER 3.0 45 dB SIZE 4 RIGHT EAR	10992557	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 45 dB SIZE 4 LEFT EAR	MINIRECEIVER 3.0 45 dB SIZE 4 LEFT EAR	10992558	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX

## ROES codes (continued)

VA ROES ITEM	DESCRIPTION	ITEM #	QTY	PRICE	COMPATIBLE PRODUCTS
MINIRECEIVER 3.0 45 DB SIZE 5 RIGHT EAR	MINIRECEIVER 3.0 45 DB SIZE 5 RIGHT EAR	2102307	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 45 DB SIZE 5 LEFT EAR	MINIRECEIVER 3.0 45 DB SIZE 5 LEFT EAR	21023211	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 60 dB SIZE 0 RIGHT EAR	MINIRECEIVER 3.0 60 dB SIZE 0 RIGHT EAR	10992559	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 60 dB SIZE 0 LEFT EAR	MINIRECEIVER 3.0 60 dB SIZE 0 LEFT EAR	10992560	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 60 dB SIZE 1 RIGHT EAR	MINIRECEIVER 3.0 60 dB SIZE 1 RIGHT EAR	10992561	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 60 dB SIZE 1 LEFT EAR	MINIRECEIVER 3.0 60 dB SIZE 1 LEFT EAR	10992562	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 60 dB SIZE 2 RIGHT EAR	MINIRECEIVER 3.0 60 dB SIZE 2 RIGHT EAR	10992563	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 60 dB SIZE 2 LEFT EAR	MINIRECEIVER 3.0 60 dB SIZE 2 LEFT EAR	10992564	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 60 dB SIZE 3 RIGHT EAR	MINIRECEIVER 3.0 60 dB SIZE 3 RIGHT EAR	10992565	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 60 dB SIZE 3 LEFT EAR	MINIRECEIVER 3.0 60 dB SIZE 3 LEFT EAR	10992566	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 60 dB SIZE 4 RIGHT EAR	MINIRECEIVER 3.0 60 dB SIZE 4 RIGHT EAR	10992567	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 60 dB SIZE 4 LEFT EAR	MINIRECEIVER 3.0 60 dB SIZE 4 LEFT EAR	10992568	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX

## ROES codes (continued)

VA ROES ITEM	DESCRIPTION	ITEM #	QTY	PRICE	COMPATIBLE PRODUCTS
MINIRECEIVER 3.0 60 dB SIZE 5 RIGHT EAR	MINIRECEIVER 3.0 60 dB SIZE 5 RIGHT EAR	21023215	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 60 dB SIZE 5 LEFT EAR	MINIRECEIVER 3.0 60 dB SIZE 5 LEFT EAR	2102316	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 65 dB SIZE 0 RIGHT EAR	MINIRECEIVER 3.0 65 dB SIZE 0 RIGHT EAR	10992569	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 65 dB SIZE 0 LEFT EAR	MINIRECEIVER 3.0 65 dB SIZE 0 LEFT EAR	10992570	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 65 dB SIZE 1 RIGHT EAR	MINIRECEIVER 3.0 65 dB SIZE 1 RIGHT EAR	10992571	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 65 dB SIZE 1 LEFT EAR	MINIRECEIVER 3.0 65 dB SIZE 1 LEFT EAR	10992572	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 65 dB SIZE 2 RIGHT EAR	MINIRECEIVER 3.0 65 dB SIZE 2 RIGHT EAR	10992573	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 65 dB SIZE 2 LEFT EAR	MINIRECEIVER 3.0 65 dB SIZE 2 LEFT EAR	10992574	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 65 dB SIZE 3 RIGHT EAR	MINIRECEIVER 3.0 65 dB SIZE 3 RIGHT EAR	10992575	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 65 dB SIZE 3 LEFT EAR	MINIRECEIVER 3.0 65 dB SIZE 3 LEFT EAR	10992576	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 65 dB SIZE 4 RIGHT EAR	MINIRECEIVER 3.0 65 dB SIZE 4 RIGHT EAR	10992577	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
MINIRECEIVER 3.0 65 dB SIZE 4 LEFT EAR	MINIRECEIVER 3.0 65 dB SIZE 4 LEFT EAR	10992578	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS PURE 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX

## ROES codes (continued)

VA ROES ITEM	DESCRIPTION	ITEM #	QTY	PRICE	COMPATIBLE PRODUCTS
MINIRECEIVER 3.0 65 dB SIZE 5 RIGHT EAR	MINIRECEIVER 3.0 65 dB SIZE 5 RIGHT EAR	21023217	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG X RIC, PURE 312T X RIC, STYLELTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
MINIRECEIVER 3.0 65 dB SIZE 5 LEFT EAR	MINIRECEIVER 3.0 65 dB SIZE 5 LEFT EAR	21023218	1	\$51.36	PURE CNG IX, CROS CNG IX, PURE CNG X RIC, PURE 312T X RIC, STYLELTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
MINIRECEIVER COMPLETE SET	MINI RECEIVER 2.0 KIT	10960146	1	\$546.00	PURE CNG IX, CROS CNG IX, PURE CNG X RIC, PURE 312T X RIC, STYLELTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
S-RECEIVER 2.0 SIZE 0 L	MINI RECEIVER 2.0 45 dB SIZE 0 LEFT EAR	10669907	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLELTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
S-RECEIVER 2.0 SIZE 0 R	MINI RECEIVER 2.0 45 dB SIZE 0 RIGHT EAR	10669908	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLELTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
S-RECEIVER 2.0 SIZE 1 L	MINI RECEIVER 2.0 45 dB SIZE 1 LEFT EAR	10669909	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLELTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
S-RECEIVER 2.0 SIZE 1 R	MINI RECEIVER 2.0 45 dB SIZE 1 RIGHT EAR	10669910	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLELTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
S-RECEIVER 2.0 SIZE 2 L	MINI RECEIVER 2.0 45 dB SIZE 2 LEFT EAR	10669911	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLELTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
S-RECEIVER 2.0 SIZE 2 R	MINI RECEIVER 2.0 45 dB SIZE 2 RIGHT EAR	10669912	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLELTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
S-RECEIVER 2.0 SIZE 3 L	MINI RECEIVER 2.0 45 dB SIZE 3 LEFT EAR	10669913	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLELTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
S-RECEIVER 2.0 SIZE 3 R	MINI RECEIVER 2.0 45 dB SIZE 3 RIGHT EAR	10669914	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLELTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
M-RECEIVER 2.0 SIZE 0 L	MINI RECEIVER 2.0 60 dB SIZE 0 LEFT EAR	10824820	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLELTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC



## ROES codes (continued)

VA ROES ITEM	DESCRIPTION	ITEM #	QTY	PRICE	COMPATIBLE PRODUCTS
M-RECEIVER 2.0 SIZE 0 R	MINI RECEIVER 2.0 60 dB SIZE 0 RIGHT EAR	10824821	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
M-RECEIVER 2.0 SIZE 1 L	MINI RECEIVER 2.0 60 dB SIZE 1 LEFT EAR	10824822	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
M-RECEIVER 2.0 SIZE 1 R	MINI RECEIVER 2.0 60 dB SIZE 1 RIGHT EAR	10824823	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
M-RECEIVER 2.0 SIZE 2 L	MINI RECEIVER 2.0 60 dB SIZE 2 LEFT EAR	10824824	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
M-RECEIVER 2.0 SIZE 2 R	MINI RECEIVER 2.0 60 dB SIZE 2 RIGHT EAR	10824825	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
M-RECEIVER 2.0 SIZE 3 L	MINI RECEIVER 2.0 60 dB SIZE 3 LEFT EAR	10824826	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
M-RECEIVER 2.0 SIZE 3 R	MINI RECEIVER 2.0 60dB SIZE 3 RIGHT EAR	10824827	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
M-RECEIVER 2.0 SIZE 4 L	MINI RECEIVER 2.0 60 dB SIZE 1 LEFT EAR	10933054	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
M-RECEIVER 2.0 SIZE 4 R	MINI RECEIVER 2.0 60 dB SIZE 1 RIGHT EAR	10933053	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
P-RECEIVER 2.0 SIZE 0 L	MINI RECEIVER 2.0 70 dB SIZE 0 LEFT EAR	10669923	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
P-RECEIVER 2.0 SIZE 0 R	MINI RECEIVER 2.0 70 dB SIZE 0 RIGHT EAR	10669924	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
P-RECEIVER 2.0 SIZE 1 L	MINI RECEIVER 2.0 70 dB SIZE 1 LEFT EAR	10669925	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC

## ROES codes (continued)

VA ROES ITEM	DESCRIPTION	ITEM #	QTY	PRICE	COMPATIBLE PRODUCTS
P-RECEIVER 2.0 SIZE 1 R	MINI RECEIVER 2.0 70 dB SIZE 1 RIGHT EAR	10669926	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
P-RECEIVER 2.0 SIZE 2 L	MINI RECEIVER 2.0 70 dB SIZE 2 LEFT EAR	10669927	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
P-RECEIVER 2.0 SIZE 2 R	MINI RECEIVER 2.0 70 dB SIZE 2 RIGHT EAR	10669928	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
P-RECEIVER 2.0 SIZE 3 L	MINI RECEIVER 2.0 70 dB SIZE 3 LEFT EAR	10669929	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
P-RECEIVER 2.0 SIZE 3 R	MINI RECEIVER 2.0 70 dB SIZE 3 RIGHT EAR	10669930	1	\$51.36	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
<b>THIN TUBES</b>					
3.0 THIN TUBE KIT	THIN TUBE 3.0 FITTING KIT	10992763	1	\$93.60	MOTION CNG X BTE , MOTION CNG P X BTE
BINAURAL THINTUBE 3.0 S0	THIN TUBE 3.0 SIZE 0	10991609	2	\$8.53	MOTION CNG X BTE , MOTION CNG P X BTE
BINAURAL THINTUBE 3.0 S1	THIN TUBE 3.0 SIZE 1	10991610	2	\$8.53	MOTION CNG X BTE , MOTION CNG P X BTE
BINAURAL THINTUBE 3.0 S2	THIN TUBE 3.0 SIZE 2	10991611	2	\$8.53	MOTION CNG X BTE , MOTION CNG P X BTE
BINAURAL THINTUBE 3.0 S3	THIN TUBE 3.0 SIZE 3	10991612	2	\$8.53	MOTION CNG X BTE , MOTION CNG P X BTE
BINAURAL THINTUBE 3.0 S4	THIN TUBE 3.0 SIZE 4	10991613	2	\$8.53	MOTION CNG X BTE , MOTION CNG P X BTE
BINAURAL THINTUBE 3.0 S0 POWER	POWER THIN TUBE 3.0 SIZE 0	10992492	2	\$8.53	MOTION CNG X BTE , MOTION CNG P X BTE
BINAURAL THINTUBE 3.0 S1 POWER	POWER THIN TUBE 3.0 SIZE 1	10992493	2	\$8.53	MOTION CNG X BTE , MOTION CNG P X BTE
BINAURAL THINTUBE 3.0 S2 POWER	POWER THIN TUBE 3.0 SIZE 2	10992494	2	\$8.53	MOTION CNG X BTE , MOTION CNG P X BTE

## ROES codes (continued)

VA ROES ITEM	DESCRIPTION	ITEM #	QTY	PRICE	COMPATIBLE PRODUCTS
BINAURAL THINTUBE 3.0 S3 POWER	POWER THIN TUBE 3.0 SIZE 3	10992495	2	\$8.53	MOTION CNG X BTE , MOTION CNG P X BTE
BINAURAL THINTUBE 3.0 S4 POWER	POWER THIN TUBE 3.0 SIZE 4	10992496	2	\$8.53	MOTION CNG X BTE , MOTION CNG P X BTE
TUBE THIN 2.0 FITTING SET	THIN TUBE 2.0 FITTING KIT	10828299	1	\$169.71	MOTION CNG NX BTE, MOTION 13 NX BTE, MOTION 13P NX BTE
THIN TUBE 2.0 1L	THIN TUBE 2.0 SIZE 1 LEFT	10828291	5	\$8.53	MOTION CNG NX BTE, MOTION 13 NX BTE, MOTION 13P NX BTE
THIN TUBE 2.0 1R	THIN TUBE 2.0 SIZE 1 RIGHT	10828292	5	\$8.53	MOTION CNG NX BTE, MOTION 13 NX BTE, MOTION 13P NX BTE
THIN TUBE 2.0 2L	THIN TUBE 2.0 SIZE 2 LEFT	10828293	5	\$8.53	MOTION CNG NX BTE, MOTION 13 NX BTE, MOTION 13P NX BTE
THIN TUBE 2.0 2R	THIN TUBE 2.0 SIZE 2 RIGHT	10828294	5	\$8.53	MOTION CNG NX BTE, MOTION 13 NX BTE, MOTION 13P NX BTE
THIN TUBE 2.0 3L	THIN TUBE 2.0 SIZE 3 LEFT	10828295	5	\$8.53	MOTION CNG NX BTE, MOTION 13 NX BTE, MOTION 13P NX BTE
THIN TUBE 2.0 3R	THIN TUBE 2.0 SIZE 3 RIGHT	10828296	5	\$8.53	MOTION CNG NX BTE, MOTION 13 NX BTE, MOTION 13P NX BTE
THIN TUBE 2.0 4L	THIN TUBE 2.0 SIZE 4 LEFT	10828297	5	\$8.53	MOTION CNG NX BTE, MOTION 13 NX BTE, MOTION 13P NX BTE
THIN TUBE 2.0 4R	THIN TUBE 2.0 SIZE 4 RIGHT	10828298	5	\$8.53	MOTION CNG NX BTE, MOTION 13 NX BTE, MOTION 13P NX BTE
<b>DOMES AND SLEEVES</b>					
EARTIP 3.0 5MM OPEN	EARTIP 3.0 S OPEN	10993678	6	\$8.53	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, MOTION CNG X BTE, MOTION CNG P BTE, Pure 312 AX, CROS 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
EARTIP 3.0 7MM OPEN	EARTIP 3.0 M OPEN	10993680	6	\$8.53	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, MOTION CNG X BTE, MOTION CNG P BTE, Pure 312 AX, CROS 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
EARTIP 3.0 10MM OPEN	EARTIP 3.0 L OPEN	10993681	6	\$8.53	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, MOTION CNG X BTE, MOTION CNG P BTE, Pure 312 AX, CROS 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX

## ROES codes (continued)

VA ROES ITEM	DESCRIPTION	ITEM #	QTY	PRICE	COMPATIBLE PRODUCTS
EARTIP 3.0 8MM TULIP	EARTIP 3.0 M TULIP	10993682	6	\$8.53	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, MOTION CNG X BTE, MOTION CNG P BTE, Pure 312 AX, CROS 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
EARTIP 3.0 12MM TULIP	EARTIP 3.0 L TULIP	10993683	6	\$8.53	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, MOTION CNG X BTE, MOTION CNG P BTE, Pure 312 AX, CROS 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
EARTIP 3.0 S ROUND DOME 1	EARTIP 3.0 S ROUND DOME 1	21022728	6	\$8.53	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, MOTION CNG X BTE, MOTION CNG P BTE, Pure 312 AX, CROS 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
EARTIP 3.0 M ROUND DOME 1	EARTIP 3.0 M ROUND DOME 1	21022729	6	\$8.53	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, MOTION CNG X BTE, MOTION CNG P BTE, Pure 312 AX, CROS 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
EARTIP 3.0 L ROUND DOME 1	EARTIP 3.0 L ROUND DOME 1	21022730	6	\$8.53	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, MOTION CNG X BTE, MOTION CNG P BTE, Pure 312 AX, CROS 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
SLEEVE 3.0 VENTED XS	SLEEVE 3.0 XS VENTED	10992863	6	\$8.53	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, MOTION CNG X BTE, MOTION CNG P BTE, Pure 312 AX, CROS 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
SLEEVE 3.0 VENTED S	SLEEVE 3.0 VENTED S TBLK	10992864	6	\$8.53	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, MOTION CNG X BTE, MOTION CNG P BTE, Pure 312 AX, CROS 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
SLEEVE 3.0 VENTED M	SLEEVE 3.0 VENTED M TBLK	10992865	6	\$8.53	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, MOTION CNG X BTE, MOTION CNG P BTE, Pure 312 AX, CROS 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
SLEEVE 3.0 VENTED L	SLEEVE 3.0 VENTED L TBLK	10992866	6	\$8.53	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, MOTION CNG X BTE, MOTION CNG P BTE, Pure 312 AX, CROS 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX

## ROES codes (continued)

VA ROES ITEM	DESCRIPTION	ITEM #	QTY	PRICE	COMPATIBLE PRODUCTS
SLEEVE 3.0 CLOSED XS	SLEEVE 3.0 CLOSED XS TBLK	10992867	6	\$8.53	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, MOTION CNG X BTE, MOTION CNG P BTE, Pure 312 AX, CROS 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
SLEEVE 3.0 CLOSED S	SLEEVE 3.0 CLOSED S TBLK	10992868	6	\$8.53	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, MOTION CNG X BTE, MOTION CNG P BTE, Pure 312 AX, CROS 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
SLEEVE 3.0 CLOSED M	SLEEVE 3.0 CLOSED M TBLK	10992869	6	\$8.53	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, MOTION CNG X BTE, MOTION CNG P BTE, Pure 312 AX, CROS 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
SLEEVE 3.0 CLOSED L	SLEEVE 3.0 CLOSED L TBLK	10992870	6	\$8.53	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, MOTION CNG X BTE, MOTION CNG P BTE, Pure 312 AX, CROS 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
SLEEVE 3.0 POWER XS	SLEEVE 3.0 POWER XS CXX	10997192	6	\$8.53	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, MOTION CNG X BTE, MOTION CNG P BTE, Pure 312 AX, CROS 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
SLEEVE 3.0 POWER S	SLEEVE 3.0 POWER S CXX	10997193	6	\$8.53	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, MOTION CNG X BTE, MOTION CNG P BTE, Pure 312 AX, CROS 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
SLEEVE 3.0 POWER M	SLEEVE 3.0 POWER M CXX	10997194	6	\$8.53	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, MOTION CNG X BTE, MOTION CNG P BTE, Pure 312 AX, CROS 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
SLEEVE 3.0 POWER L	SLEEVE 3.0 POWER L CXX	10997195	6	\$8.53	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, MOTION CNG X BTE, MOTION CNG P BTE, Pure 312 AX, CROS 312 AX, STYLETTO IX/AX, CROS STYLETTO IX/AX
CLICK SLEEVES XS VENTED	CLICK SLEEVE 2.0 XS VENTED	10951152	6	\$8.53	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC

## ROES codes (continued)

VA ROES ITEM	DESCRIPTION	ITEM #	QTY	PRICE	COMPATIBLE PRODUCTS
CLICK SLEEVES S VENTED	CLICK SLEEVE 2.0 S VENTED	10951146	6	\$8.53	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
CLICK SLEEVES M VENTED	CLICK SLEEVE 2.0 M VENTED	10951148	6	\$8.53	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
CLICK SLEEVES L VENTED	CLICK SLEEVE 2.0 L VENTED	10951150	6	\$8.53	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
CLICK SLEEVES S CLOSED	CLICK SLEEVE 2.0 S CLOSED	10951147	6	\$8.53	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
CLICK SLEEVES M CLOSED	CLICK SLEEVE 2.0 M CLOSED	10951149	6	\$8.53	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
CLICK SLEEVES L CLOSED	CLICK SLEEVE 2.0 L CLOSED	10951151	6	\$8.53	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
CLICK DOMES OPEN 4MM	CLICK DOME S OPEN	10426016	6	\$5.12	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
CLICK DOMES OPEN 8MM	CLICK DOME M OPEN	10426014	6	\$5.12	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
CLICK DOMES OPEN 10MM	CLICK DOME L OPEN	10426015	6	\$5.12	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
CLICK DOMES CLOSED 6MM	CLICK DOME S CLOSED	10426017	6	\$5.12	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
CLICK DOMES CLOSED 8MM	CLICK DOME M CLOSED	10426012	6	\$5.12	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
CLICK DOMES CLOSED 10MM	CLICK DOME L CLOSED	10426013	6	\$5.12	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC

## ROES codes (continued)

VA ROES ITEM	DESCRIPTION	ITEM #	QTY	PRICE	COMPATIBLE PRODUCTS
CLICK DOMES DOUBLE 8/10MM	CLICK DOME S DOUBLE	10426018	6	\$5.12	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
CLICK DOMES DOUBLE 10/12MM	CLICK DOME M DOUBLE	10426019	6	\$5.12	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
DOME 8/12MM SEMI OPEN CLICK RIC	CLICK DOME M SEMI OPEN	10670790	6	\$5.12	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
LIFE TIP 8 MM OPEN	LIFE TIP REFILL 8MM	10054916	1	\$8.53	MOTION CNG NX BTE, MOTION 13 NX BTE. MOTION 13P NX BTE
LIFE TIP 8 MM CLOSED	LIFE TIP REFILL 8MM FULL	10169330	1	\$8.53	MOTION CNG NX BTE, MOTION 13 NX BTE. MOTION 13P NX BTE
LIFE TIP 4 MM	LIFE TIP REFILL 4 MM	10169328	1	\$8.53	MOTION CNG NX BTE, MOTION 13 NX BTE. MOTION 13P NX BTE
LIFE TIP 6 MM	LIFE TIP REFILL 6MM	10169329	1	\$8.53	MOTION CNG NX BTE, MOTION 13 NX BTE. MOTION 13P NX BTE
LIFE TIP 10 MM	LIFE TIP REFILL 10MM	10054917	1	\$8.53	MOTION CNG NX BTE, MOTION 13 NX BTE. MOTION 13P NX BTE
DOUBLE TIP 8/10	LIFE TIP DOUBLE REFILL 8/10MM	10175566	1	\$8.53	MOTION CNG NX BTE, MOTION 13 NX BTE. MOTION 13P NX BTE
DOUBLE TIP 10/12	LIFE TIP DOUBLE REFILL 10/12MM	10175715	1	\$8.53	MOTION CNG NX BTE, MOTION 13 NX BTE. MOTION 13P NX BTE
<b>CHARGERS</b>					
QI CHARGING PAD		10245297			
CLINIC MULTI-CHARGER		21025145	1	\$79	ALL HEARING AIDS THAT HAVE GALVANIC CHARGING
CLINIC MULTI-CHARGER STAND		21027726	1	\$15	ALL HEARING AIDS THAT HAVE GALVANIC CHARGING
INSIO CNG AX CHARGER	CHARGER FOR INSIO CUSTOM HEARING AIDS	21025145	1	\$175.76	INSIO CNG AX

## ROES codes (continued)

VA ROES ITEM	DESCRIPTION	ITEM #	QTY	PRICE	COMPATIBLE PRODUCTS
PURE CNG PORTABLE CHARGER	PORTABLE CHARGER FOR PURE C&G IX/AX	10993118	1	\$175.76	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX
PURE CNG DRY & CLEAN CHARGER	DRY&CLEAN CHARGER FOR PURE C&G IX/AX	10991650	1	\$218.40	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX
PURE CNG STANDARD CHARGER	CHARGER FOR PURE C&G IX/AX	10996345	1	\$175.76	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX
CHARGER STYLETTO IX	CHARGER FOR STYLETTO IX	10997503	1	\$175.76	STYLETTO IX, CROS STYLETTO IX
CHARGER STYLETTO AX	CHARGER FOR STYLETTO AX	21003216	1	\$175.76	STYLETTO AX RIC, CROS STYLETTO AX
DRY&CLEAN CHARGER STYLETTO IX/AX	DRY&CLEAN CHARGER FOR STYLETTO IX/AX	21007593	1	\$218.40	STYLETTO IX/AX, CROS STYLETTO IX/AX
DRY & CLEAN CHARGER MOTION X BTE P	DRY & CLEAN CHARGER FOR MOTION CNG P X BTE	10991638	1	\$218.40	MOTION CNG X P BTE
DRY & CLEAN CHARGER MOTION X BTE SP	DRY & CLEAN CHARGER FOR MOTION CNG SP X BTE	10991644	1	\$218.40	MOTION CNG X SP BTE
CHARGER MOTION XBTE P	CHARGER FOR MOTION CNG P X BTE	10996335	1	\$175.76	MOTION CNG P BTE
CHARGER MOTION X BTE SP	CHARGER FOR MOTION CNG SP X BTE	10996340	1	\$175.76	MOTION CNG SP BTE
INDUCTIVE CHARGER II	INDUCTIVE CHARGER II FOR PURE AND MOTION	10966942	1	\$165.57	PURE CNG X RIC, MOTION CNG X, MOTION CNG NX
CHARGER STYLETTO X	CHARGER FOR STYLETTO X	10972015	1	\$175.76	STYLETTO X RIC
CHARGER ACTIVE PRO	CHARGER FOR ACTIVE PRO	10991347	1	\$175.76	ACTIVE PRO
CHARGER INDUCTIVE SET LI-ION	LI-ION INDUCTIVE CHARGER SET	10948199	1	\$175.76	CELLION PX, PURE CNG NX, MOTION CNG NX



## ROES codes (continued)

VA ROES ITEM	DESCRIPTION	ITEM #	QTY	PRICE	COMPATIBLE PRODUCTS
<b>WAXGUARDS</b>					
WAXGUARD 3.0 NANOCARE	EARWEAR 3.0 NANOCARE WAXGUARD	10993607	1	\$0.57	PURE CNG AX, INSIO AX, MOTION CNG X BTE, MOTION CNG P X BTE, STYLETTO AX
CERUSTOP XL	WAXGUARD FOR HP EAR-MOLD	10999334	1	\$0.00	PURE CNG AX HP MOLD
WAXGUARD RIC MINI	MINI RECEIVER WAXGUARD	10966257	10	\$2.98	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
QUICKGUARD PACK	CXX QuickGuard	10962504	1	\$3.41	INSIO NX CUSTOMS
WAXGUARD C-GRID	C-GRID WAXGUARD	10423507	16	\$8.53	INSIO NX CUSTOMS
WAX SPRING,RCVR TUBE ITE	WAX SPRING FOR RCVR TUBE	07080208	1	\$0.00	INSIO NX CUSTOMS
<b>CONCHA LOCKS</b>					
CONCHA LOCK 3.0 SLEEVES S RECEIVERS		21024885	10	\$0.00	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS 312 AX, STYLETTO AX, CROS STYLETTO AX, MOTION CNG X BTEs
CONCHA LOCK 3.0 SLEEVES M RECEIVERS		21024886	10	\$0.00	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS 312 AX, STYLETTO AX, CROS STYLETTO AX, MOTION CNG X BTEs
CONCHA LOCK 3.0 SLEEVES P RECEIVERS		21024887	10	\$0.00	PURE CNG IX, CROS CNG IX, PURE CNG AX, CROS CNG AX, PURE 312 AX, CROS 312 AX, STYLETTO AX, CROS STYLETTO AX, MOTION CNG X BTEs
RECEIVER CONCHA LOCK S45dB	MINI RECEIVER 2.0 CONCHA LOCK FOR S RECEIVER	10601877	10	\$0.00	PURE CNG X RIC, PURE 312T X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
RECEIVER CONCHA LOCK M60dB	MINI RECEIVER 2.0 CONCHA LOCK FOR M RECEIVER	10937826	10	\$0.00	PURE CNG X RIC, PURE 312T X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC

## ROES codes (continued)

VA ROES ITEM	DESCRIPTION	ITEM #	QTY	PRICE	COMPATIBLE PRODUCTS
RECEIVER CONCHA LOCK P70db	MINI RECEIVER 2.0 CONCHA LOCK FOR P RECEIVER	10673516	10	\$0.00	PURE CNG X RIC, PURE 312T X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
<b>EARHOOKS</b>					
EARHOOK MOTION C-G NX	EARHOOK FOR MOTION CNG NX BTE	10959864	1	\$0.00	MOTION CNG NX BTE
MOTION PX/NX EARHOOK	EARHOOK FOR MOTION 13 NX BTE	10947168	1	\$0.00	MOTION 13 NX BTE
THINTUBE ADAPTER MOTION	THIN TUBE ADAPTER FOR MOTION CNG NX BTE	10958343	1	\$3.12	MOTION CNG NX BTE
MOTION PX/NX TUBE ADAPTER	THIN TUBE ADAPTER FOR MOTION 13 NX BTE	10947303	1	\$0.00	MOTION 13 NX BTE
EARHOOK MOTION X SMALL	EARHOOK MOTION X SMALL FOR	10980403	1	\$0.00	MOTION CNG X BTE, MOTION CNG P X BTE
EARHOOK MOTION X STANDARD	EARHOOK MOTION X STANDARD	10980402	1	\$0.00	MOTION CNG X BTE, MOTION CNG P X BTE
EARHOOK MOTION X SP SMALL	EARHOOK MOTION X SP SMALL	10993620	1	\$0.00	MOTION CNG SP X BTE
<b>MIC COVERS</b>					
CUSTOM MIC COVER BEIGE, TAN, MOCHA	MIC COVERS FOR CUSTOM HEARING AIDS - BEIGE, TAN, MOCHA	10947542	1	\$0.00	INSIO C&G AX, INSIO NX
CUSTOM MIC COVER BROWN	MIC COVER FOR CUSTOM HEARING AIDS - BROWN	10947776	1	\$0.00	INSIO C&G AX, INSIO NX

## ROES codes (continued)

VA ROES ITEM	DESCRIPTION	ITEM #	QTY	PRICE	COMPATIBLE PRODUCTS
CUSTOM MIC COVER DARK BRN	MIC COVER FOR CUSTOM HEARING AIDS - DARK BROWN	21016202	1	\$0.00	INSIO C&G AX, INSIO NX
CUSTOM MIC COVER BLK	MIC COVER FOR CUSTOM HEARING AIDS - BLACK	21016200	1	\$0.00	INSIO C&G AX, INSIO NX
<b>MISC TOOLS</b>					
3.0 MEASURING TOOL	EARWEAR 3.0 MEASURING TOOL	10991403	1	\$0.00	PURE CNG AX, CROS CNG AX, MOTION CNG X BTE, MOTION CNG P X BTE
EAR-TIP 3.0 INSERTION TOOL	EARWEAR 3.0 INSERTION TOOL	10994214	1	\$0.00	PURE CNG AX, CROS CNG AX, MOTION CNG X BTE, MOTION CNG P X BTE
TOOL MEASURING	RECEIVER MEASURING TOOL	10174571	1	\$0.00	PURE CNG X RIC , PURE 312T X RIC, STYLETTO X RIC. PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
MINI RECEIVER TOOL	MINI RECEIVER REMOVAL TOOL	10674370	1	\$0.00	PURE CNG X , PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX, PURE 312T NX RIC, PURE 13 NX RIC
THIN TUBE 3.0 PIN EJECTOR TOOL	EARHOOK REMOVAL TOOL	10600047	1	\$0.00	MOTION CNG X BTE, MOTION CNG P X BTE
TOOL CLICK MOLD	CLICK MOLD SERVICE TOOL	10597930	1	\$0.00	PURE CNG X RIC, PURE 312T X RIC, STYLETTO X RIC, PURE CNG NX RIC, PURE 312T NX RIC, PURE 13 NX RIC
EARTIP 3.0 HANDLING TOOL	3.0 EARMOLD SERVICE TOOL	10984033	1	\$0.00	MOTION CNG X BTE, MOTION CNG P X BTE
<b>ACCESSORY PARTS &amp; STREAMLINE TV AND MIC</b>					
STREAMLINE TV POWER SUPPLY	POWER SUPPLY FOR STREAM-LINE TV	10944815	1	\$0.00	STREAMLINE TV
SMARTCONNECT POWER ADAPT SET	SMARTCONNECT POWER ADAPT SET	10824346	1	\$0.00	SMARTCONNECT



**Our mission is to help improve the quality of life of Veterans,  
Active Military and American Indians with better hearing.**

**Signia Government Services: (800) 955-4327 | Consumer Hotline: (800) 350-6093**

The Bluetooth® word mark and logos are owned by the Bluetooth SIG, Inc. and any use of such marks by Sivantos is under license. Other trademarks and trade names are those of their respective owners.

Android is a trademark of Google Inc.

iOS is a trademark or registered trademark of Cisco in the U.S. and other countries and is used under license.

"Made for iPhone" means that an electronic accessory has been designed to connect specifically to iPhone and has been certified by the developer to meet Apple performance standards. Apple is not responsible for the operation of this device or its compliance with safety and regulatory standards. Please note that the use of this accessory with iPhone may affect wireless performance.



Signia is a registered trademark of WSAUD A/S. 4/24 SIV-18645

**Be  
Brilliant™**